

**TABLE 6.2D: RESUSCITAIRE® RADIANT WARMER, 120V, 240V VHA CART
ASSEMBLY, VARITILT BASSINET, PARTS LIST**

(Sheet 1 of 3)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
1		OUTER LOWER COLUMN	82 011 00
2		SWIVEL CASTER 5.00, W/BRAKE, GRAY	81 005 07
3		BASE CASTING	81 020 01
4		BASSINET CRADLE WELDMENT	82 011 01
5		FRONT HANDLE	82 011 05
6		PIVOT SUPPORT	82 001 07
7		NOT USED	
8		LATCH BLOCK	81 020 17
9		PASS-THRU DRAWER	81 020 08
10		RESUS COMPARTMENT REAR PLATE	81 020 24
11		PASS-THRU DRAWER SLIDE	81 101 12
12		RESUS COMPARTMENT COVER	81 020 18
13		LEAD SCREW	82 011 03
14		STANDOFF, 5/8 H X 10-32 ID 1.25 L S	82 011 04
15		SUCTION BOTTLE HOLDER	81 000 03
16		RESUS COMPARTMENT HINGE	81 020 22
17		CRANK KNOB ASSEMBLY	82 011 10
18		CORD WITH NYLON EYELETS	81 020 23
19		SUCTION BOTTLE 800 CC	81 001 15
20		NYLON NIPPLE 3/16 X 3/16 TUBE	81 000 13
21		PLASTIC CLIP SNAP, SLFADH	78 163 40
22		RIGHT BASSINET SHROUD	82 011 11
23		SUCTION TUBING W/END .25ID 42.0" LG	82 010 10
24		SUCTION TUBING W/END .25ID 6' LG	81 001 16
25		BASSINET ASSEMBLY	82 011 20
26		BASE DRESS PLATE, VHA	82 011 00
27		NOT USED	
28		SUCTION FLOW LABEL	81 101 23
29		120V ACTUATOR	82 010 40
		220V ACTUATOR	82 010 41
30		SCREW, 4-40 X 1/4 TR PH SS NYLOK	99 010 67
31		SCREW, 8-32 X 3/8 FL PH SS 100D NYK	99 031 54
32		SCREW, 4-40 X 1/4, FL PH SS	99 010 62
33		SCREW, 8-32 X 3/8 TR PH SS NYLOK	99 031 52
34		SCREW, 6-32 X 1/4 PN PH S Z1 SEMSC	99 022 84

Change 5

**TABLE 6.2D: RESUSCITAIRE® RADIANT WARMER, 120V, 240V VHA CART
ASSEMBLY, VARITILT BASSINET, PARTS LIST**

(Sheet 2 of 3)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
35		SCREW, 10-32 X 1/2 TR PH SS NYLOK	99 042 05
36		SCREW, 4-40 X 3/8 RD PH SS NYLOK	99 011 14
37		SCREW, 4-40 X 5/8 CP SK BUTTON SS	99 056 66
38		HEX NUT, 3/8-16 ES S ZI	99 111 25
39		NOT USED	
40		SPLIT NUT, SPECIAL 1/2-13	78 156 20
41		CRADLE LINK	78 156 25
42		SCREW, 6-32 X 1/4 RD PH SS NYLOK	99 022 86
43		SCREW, M6 X 1.0 X 0.35 LG CP SK SS	99 998 76
44		ACORN NUT, 8-32 SS	99 106 21
45		SCREW, 3/8-16 X 4.00 CP HX S ZI GR5	99 068 70
46		LOWER INNER COLUMN	82 010 13
47		LOCK WASHER, INT, #8 SS	99 122 92
48		FLAT WASHER, #10 SS 0.062 THICK	99 123 62
49		AC POWER RECEPTACLE	17 735 63
50		NYLATCH ROLLER CATCH	81 100 63
51		LOCK WASHER, INT, 1/4 SS	99 125 54
52		LOCTITE #242	AR
53		NYLON PLUG BUTTON, 0.88 HOLE	78 161 20
54		HEX NUT, 8-32 KEPS S ZI	99 106 32
55		AC POWER RECEPTACLE, FEMALE 10A	17 735 64
56		ROCKER SWITCH, DPDT 10A, 250VAC	17 638 52
57		THERMAL CIRCUIT BREAKER, SP 6A 120V	17 BH 236
		THERMAL CIRCUIT BREAKER, SP 3A 220V	17 BH 231
58		DPDT RELAY, 10A 250VAC, UL/CSA	17 653 09
		DPDT RELAY, 10A 250VAC, VDE	17 653 10
59		SUCTION JAR REPLACEMENT LABEL	81 101 24
60		RELAY MOUNTING BRACKET	82 010 06
61		OUTER COLUMN SEPARATOR PLATE	82 010 04
62		SEPARATOR PLATE GASKET	82 010 07
63		RW VHA WIRE HARNESS	82 010 05
64		CORD LOCK ADAPTER	17 734 12
65		GROUND SYMBOL LABEL	68 212 00
66		RTV COMPOUND, CLEAR, DC732	AR
67		ACTUATOR WIRE HARNESS	82 010 11

Change 5

**TABLE 6.2D: RESUSCITAIRE® RADIANT WARMER, 120V, 240V VHA CART
ASSEMBLY, VARITILT BASSINET, PARTS LIST**

(Sheet 3 of 3)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
68		POWER IN/POWER OUT WIRE HARNESS	82 010 12
69		NOT USED	
70		NOT USED	
71		COVER MUST BE REPLACED LABEL	82 010 14
72		SCREW, M6 X 6.0 LG CP SK SS	99 998 77
73		NOT USED	
74		ACTUATOR MOUNTING SPACER	82 010 30
75		SCREW, 10-32 X 3/8 TR PH SS NYLOK	99 041 36
76		TRANSPORT HEIGHT WARNING LABEL	82 000 19
77		HEX STANDOFF, 1/4-20 OD X 1.12L	78 156 30
78		LOCK WASHER, SP, 6 MM SS	99 999 17
79		NON-METALLIC FLANGE BEARING 0.38 ID	68 121 35
80		BRONZE FLANGE BEARING 0.38 ID, 0.50 OD	78 156 65
81		NEOPRENE BUMPER, 0.25 ID, 0.62 OD, 1.0 LG	81 101 16
82		SCREW, 10-32 X 5/8 CP DSK SS	99 042 55
83		SCREW, 10-32 X 3/4 CP SK SS NYLOK	99 042 91
84		SCREW, 10-32 X 2.00 CP HX SS	99 045 24
85		NOT USED	
86		SCREW, 1/4-20 X 1.12 CP SK SS	99 058 02
87		SELF LOCKING HEX NUT, 1/2 - 13 S ZI	99 113 01
88		LOCK WASHER, EXT #8 SS	99 122 93
89		LOCK WASHER, SP 310 S CA	99 124 16
90		NOT USED	
91		SPACER, 0.257 ID, 0.38 OD, 0.62 LG SS	99 125 68
92		SPACER, 0.386 ID, 0.62 HEX, 0.38 LG SS	99 126 75
93		NOT USED	
94		LOCTITE #411	AR
95		LEFT BASSINET SHROUD	82 011 09
96		CANOE CLIP, PANEL PLUG	78 155 55
97		SCREW, 10-32 X 1/2 FL PH SS	99 042 02
98		SHOULDER SCREW, SPECIAL, 14/ -20	82 011 30

Change 5

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

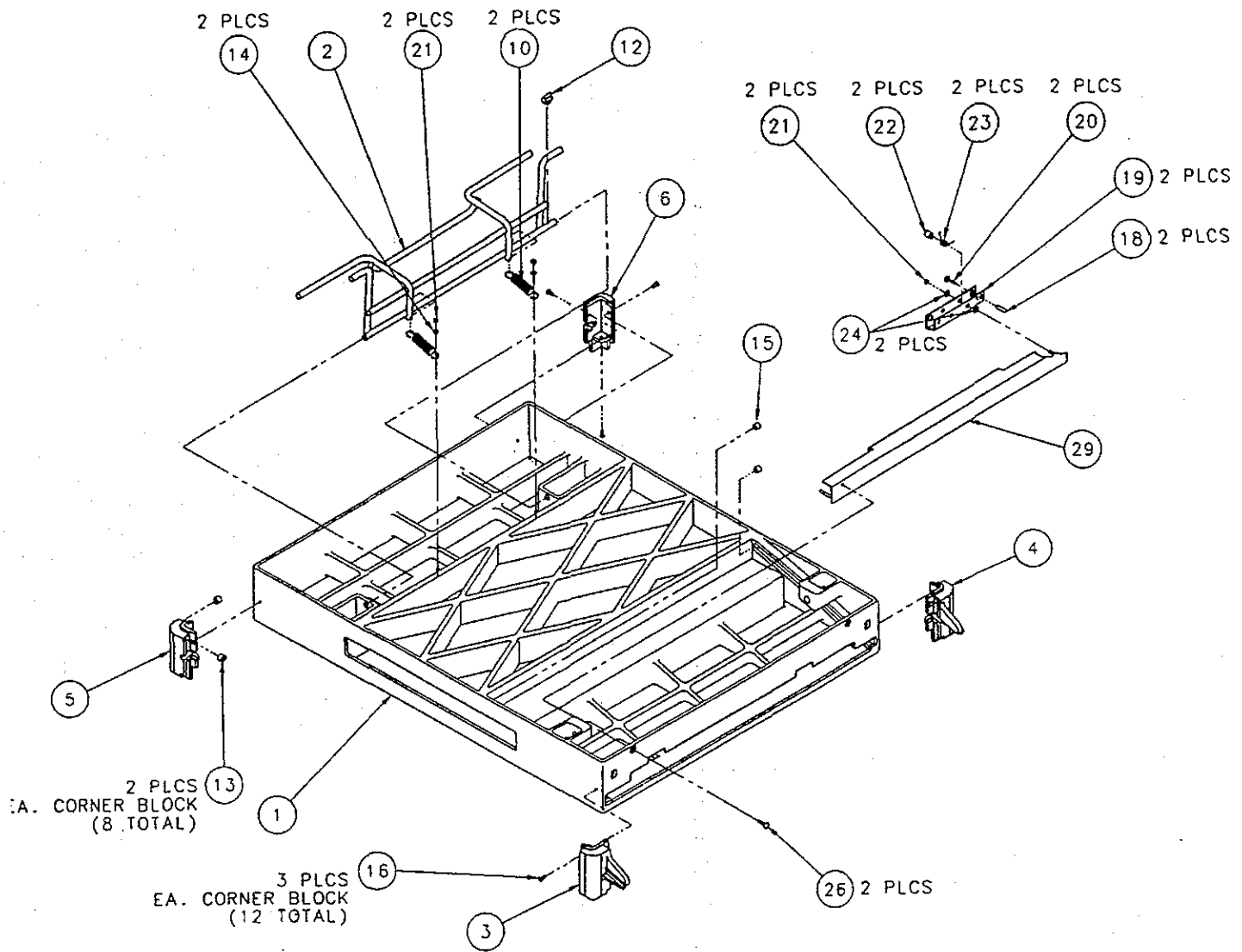


FIGURE 6.3: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
BASSINET SUB-ASSEMBLY
(SHEET 1 OF 1)

**TABLE 6.3: RESUSCITAIRE® RADIANT WARMER BASSINET SUB-ASSEMBLY,
PARTS LIST**

(Sheet 1 of 1)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		BASSINET SUB-ASSEMBLY	81 100 72
1		SUPPORT, MTRS, RESUS	81 020 07
2		HANDLE, BASSINET TILT	82 000 10
3		CORNER BLOCK, PANEL MTG, R FRONT	81 100 03
4		CORNER BLOCK, PANEL MTG, L FRONT	81 100 04
5		CORNER BLOCK, PANEL MTG, R REAR	81 100 46
6		CORNER BLOCK, PANEL MTG, L REAR	81 100 47
7THRU 9		NOT USED	
10		SPRING, EXT, .50 OD, .055W, 1.75 LG	81 100 34
11		NOT USED	
12		CLIP, SPRING, .25-.31DIA .50LG	78 162 65
13		TUBING, CLR PVC, 1/4 ID .270 LG	08 136 03
14		WSHR, FL, .19ID X .75 OD .06T S CA	99 123 32
15		TUBING, CLR PVC, 3/16ID 0.27LG	22 025 49
16		SCR, 6-32X3/8 PN PH SS NYLOK	99 023 54
17		NOT USED	
18		CAP, TUBING, VINYL, .172ID .50LG	81 100 68
19		BRACKET, AL CHNL, FRT WALL, REL	81 100 67
20		SCR, 10-32x3/4 CP SK SS NYLOK	99 042 91
21		NUT, HEX, 6-32 SS NYLOK	99 105 35
22		SPACER, .200ID .38OD .25LG CPVC	99 124 41
23		SPRING, TOR, .45W SS	81 100 69
24		WSHR, SHLDR, .14ID .22DS .05LS N	99 122 46
25		NOT USED	
26		SCR, 6-32X7/16 TR PH SS NYLOK	99 023 64
27 AND 28		NOT USED	
29		BASSINET SUPPORT	81 020 31

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

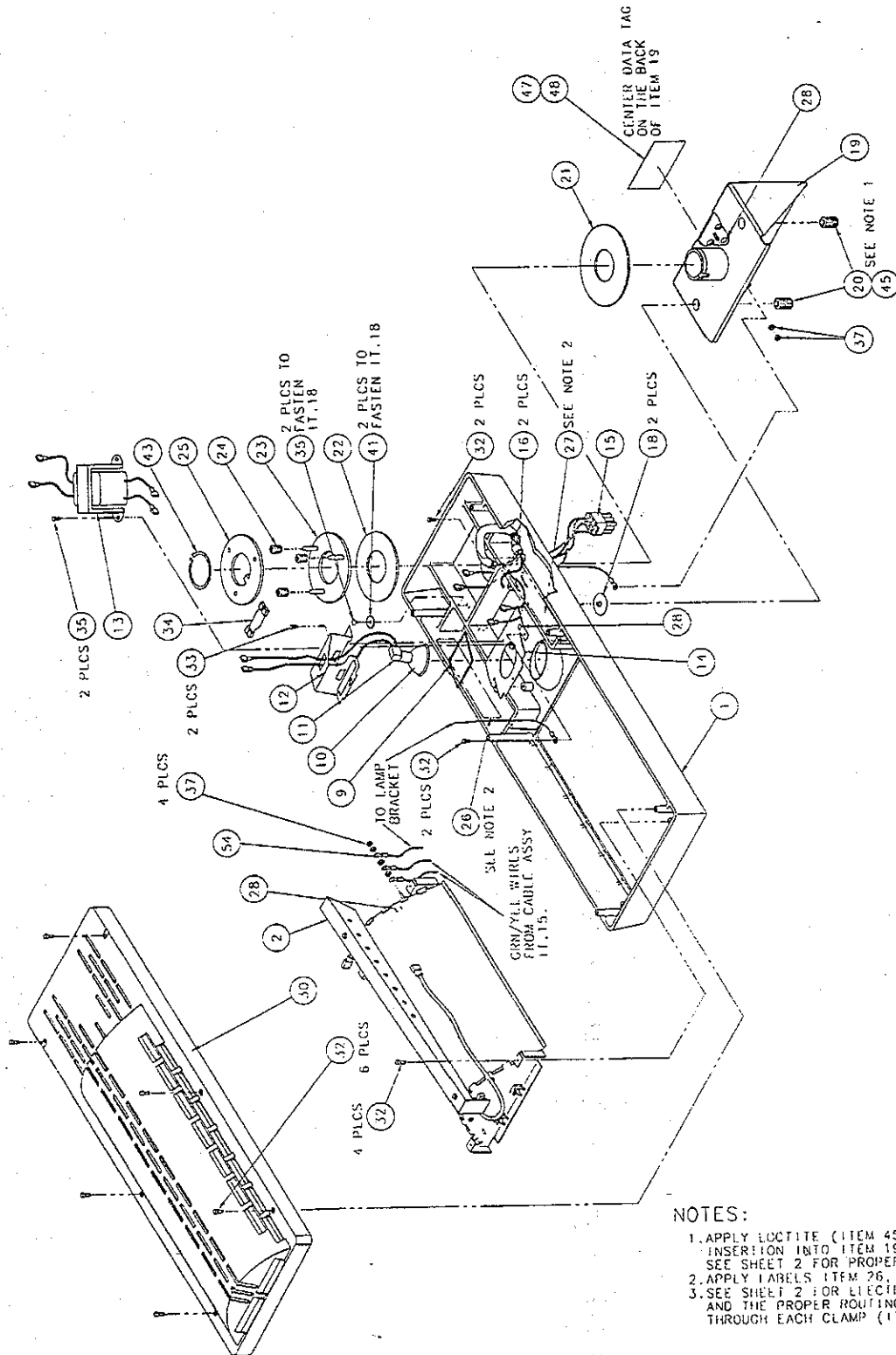


**FIGURE 6.4: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
WARMER HEAD ASSEMBLY
(SHEET 1 OF 1)
Change 1**

**TABLE 6.4: RESUSCITAIRE® RADIANT WARMER,
WARMER HEAD ASSEMBLY,
PARTS LIST**

(Sheet 1 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		WARMER HEAD ASSEMBLY, 120 V	81 200 71
		WARMER HEAD ASSEMBLY, 220/240 V	81 200 81
		WARMER HEAD ASSEMBLY, 230 V	81 200 84
1		CHASSIS, SHROUD 	
		(Order label item 11. Refer to Figure 6.1.)	81 020 10
2	Refer to Table 6.5	WARMER HD, HEATER SUB-ASSY, 120 V	81 200 75
		WARMER HD, HEATER SUB-ASSY, 220-240 V	81 200 85
3 THRU 8		NOT USED	
9		WINDOW, GLASS, SQUARE	78 260 36
10		LAMP, INCAND, QTZ-HALOGEN, 12 V 50 W	17 807 41
11		CONNECTOR ASSY, EXAM LIGHT	78 258 70
12		BRACKET, LAMP MOUNTING	78 260 25
13		TRANSFORMER ASSY, EXAM LT, 120 V	78 265 95
		TRANSFORMER ASSY, EXAM LT, 240 V	81 203 80
14		BRACKET, LAMP	81 020 28
15		CABLE ASSY, ELEC MDL TO WRMR HD	81 020 30
16		CLAMP, CABLE, LOOP TY, NYL,.375ID	17 725 64
17		NOT USED	
18		KEEPER, BALL PLUNGER	81 200 19
19		BRACKET, HEATER HEAD PIVOT	81 020 11
20		BALL PLUNGER, 5/8-11 THRD	81 200 17
21		WASHER, FL, 2.02OD 4.75OD .19PLETH	81 020 40
22		WASHER, THRUST, 2.02ID 4.25OD	81 200 15
23		WASHER, FL, BOT, W/STUDS, 4.250OD	81 200 22
24		SPRING, CPRS .480OD, .06x.09W 1" L	81 200 13
25		WASHER,FL, TOP, 2.001ID X 4.250OD	81 200 18
26		LABEL, HEATER REPL, 120 V	81 200 34
		LABEL, HEATER REPL, 220 - 240 V	81 200 35
27		LABEL, LAMP REPL EN, SP, FR, GR, IT	81 200 33
28		LABEL, GROUND SYMBOL, PROTECTIVE	68 212 00
29		NOT USED	
30		TOP COVER, WARMER 	
		(Order label item 4. Refer to Figure 6.1.)	81 020 26

Change 1

**TABLE 6.4: RESUSCITAIRE® RADIANT WARMER,
WARMER HEAD ASSEMBLY,
PARTS LIST**

(Sheet 2 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
31		NOT USED	
32		SCR, 6-32x3/8 TR PH SS NYLOK	99 023 51
33		SCR, 6-32X1/4 PN PH SS, SEMS EX	99 022 83
34		FILTER ASSY SERIES 01 ONLY	81 200 39
35		SCR, 8-32X1/4 TR PH SS NYLOK	99 030 76
36		NOT USED	
37		NUT, HEX, 6-32 KEPS S ZI	99 105 34
38		WIRE TIE, 4.5"	12 995 00
39 AND 40		NOT USED	
41		WSHR, FL, 19ID x .750D x .06T S CA	99 123 32
42		NOT USED	
43		RING, RTNG, EXT, S (WAL5100-200-ZD)	99 187 51
44		NOT USED	
45		LOCTITE #242	99 901 77
46		NOT USED	
47		LABEL, DATA TAG	81 001 01
48 THRU 53		NOT USED	
54		GROUND WIRE, #6 RING, 4.00LG	81 020 41

Change 1

This page is intentionally left blank.

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

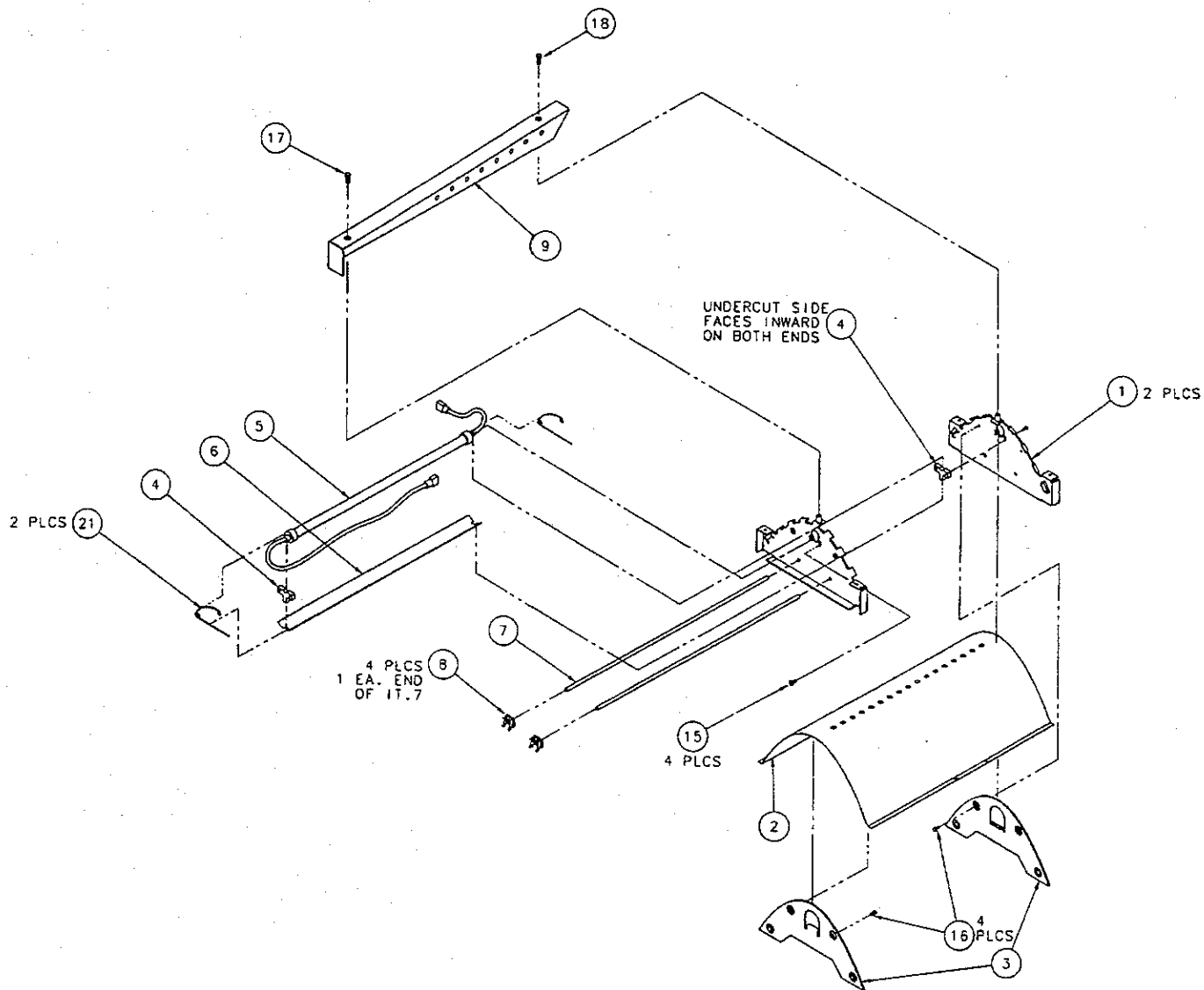


FIGURE 6.5: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
WARMER HEAD SUBASSEMBLY
(SHEET 1 OF 1)

**TABLE 6.5: RESUSCITAIRE® RADIANT WARMER,
WARMER HEAD SUBASSEMBLY,
PARTS LIST**

(Sheet 1 of 1)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		WARMER HEAD HEATER SUBASSEMBLY, 120 V	81 200 75
		WARMER HEAD HEATER SUBASSEMBLY, 220/240 V . .	81 200 85
1		END PLATE, MAIN REFLECTOR	81 020 27
2		REFLECTOR, MAIN	81 200 08
3		BRACKET, PARABOLA	81 200 40
4		BRACKET, HEATER MOUNTNG	81 200 42
5		HEATER REPLACEMENT KIT, QUARTZ, 120 V	81 990 50
		HEATER REPLACEMENT KIT, QUARTZ, 220V-240 V	81 990 51
6		REFLECTOR, SECONDARY	81 200 24
7		ROD, HEATER GUARD	81 200 07
8		CLIP, RETAINING, 3/16 SHAFT	81 200 29
9		HEAT SHIELD	81 020 29
10 THRU 14		NOT USED	
15		SCR, 6-32X1/4 PN PH SS, SEMS EX	99 022 83
16		SCR, 6-32 X 3/8 TR PH SS	99 023 31
17		SCR, 6-32X 3/8 TR PH SS NYLOK	99 023 51
18		SCR, 6-32 X1.00 TR PH SS NYLOK	99 025 23
19 AND 20		NOT USED	
21		SPRING, RETAINING	81 200 43

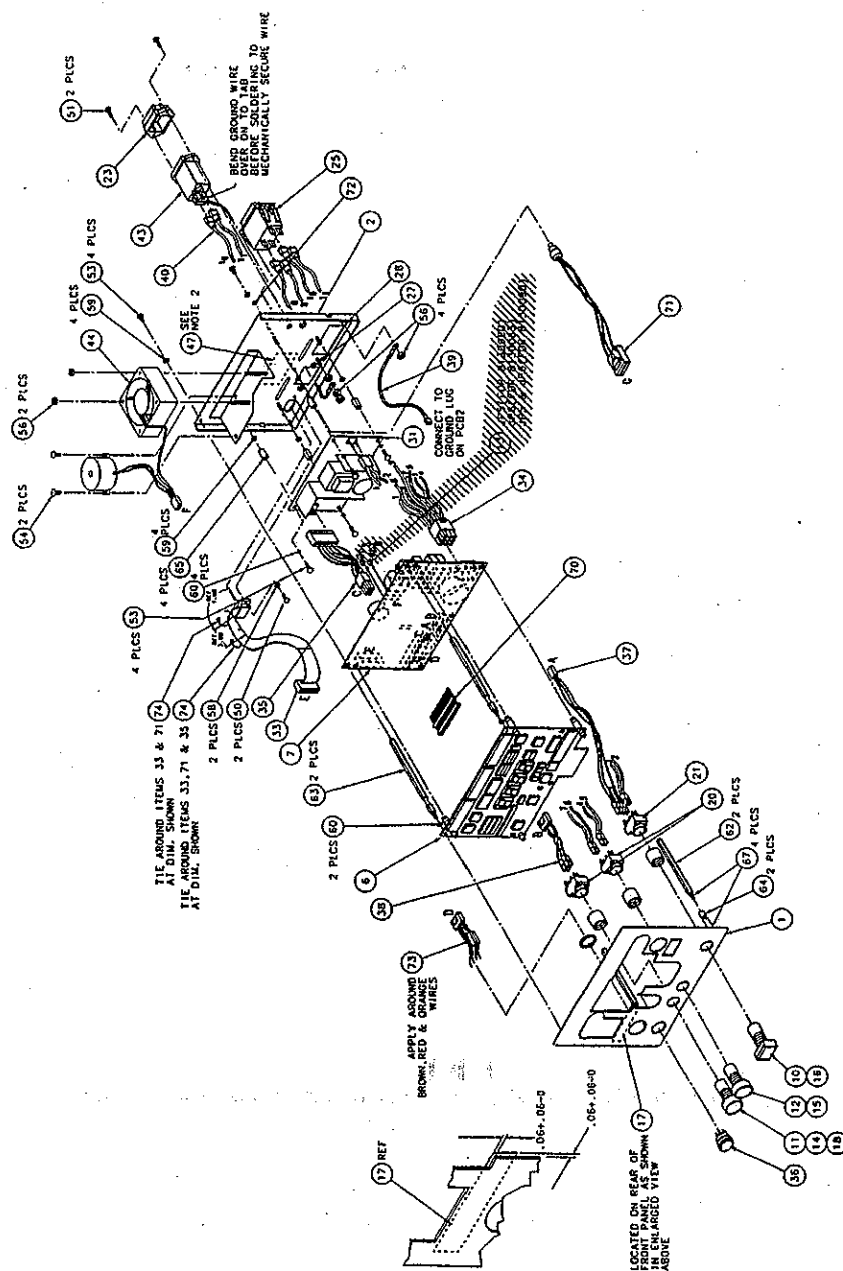


FIGURE 6.6: PARTS LOCATION DIAGRAM,
RESUSCITAIRE® RADIANT WARMER
CONTROLLER ASSEMBLY
(SHEET 1 OF 1)

**TABLE 6.6: RESUSCITAIRE® RADIANT WARMER CONTROLLER ASSEMBLY,
PARTS LIST**

(Sheet 1 of 3)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		CONTROLLER ASSEMBLY, 120V	81 300 70
		CONTROLLER ASSEMBLY, 220-240V	81 300 80
		CONTROLLER ASSEMBLY, 220V	81 300 83
1		PANEL, FRONT , ELECTRICAL MODULE 	
		(Order label item 12. Refer to Figure 6.1.)	81 300 11
2		PANEL, REAR, ELECTRICAL MODULE 	
		(Order label item 10. Refer to Figure 6.1.)	81 300 13
3 THRU 5		NOT USED	
6	Refer to Table 6.14	PCB1 ASSY, DISPLAY BOARD W/BABY MODE	X81 307 75
7	Refer to Table 6.15	PCB2 ASSY, PWR & CNTRL BOARD, 120V	X81 308 70
		PCB2 ASSY, PWR & CNTRL BOARD, 220-240V	X81 308 80
		PCB2 ASSY PWR & CONTROL 220V	X81 308 83
8 AND 9		NOT USED	
10		ACTUATOR, SW, PB, LIGHT, RECT, MNT.	17 683 01
11		ACTUATOR, SW, PB, LIGHT, RD, MOM.	17 683 08
12		ACTUATOR, SW, PB, LIGHT, RD, MNT.	17 683 09
13		NOT USED	
14		LENS, RD, FL, XPARENT, KEY INSERT	17 682 98
15		LENS, RD, FL, XPARENT, LAMP INSERT	17 682 99
16		LENS, RECT, FLAT, XLUCENT, WHITE	17 683 32
17		INSULATOR, MYLAR, 2.25 x .38x .004	81 000 08
18		LAMP, INCAND, CLR, T-1 3/4, 14V	17 683 20
19		NOT USED	
20		CONTACT BLK, SW, PB, DPST, N.O.	17 683 11
21		CONTACT BLK, SW, PB, 3PST, 2NO/NC	17 683 18
22		NOT USED	
23		CORD LOCK ADAPTER	17 734 12
24		NOT USED	
25		CKT BRKR, DPST, 12.0A 250V	17 BH2 12
		CKT BRKR, DPST, 6.0A 250V	17 BH2 09
26		NOT USED	
27		LABEL, GROUND SYMBOL, PROTECTIVE	68 212 00
28		LABEL, GROUND SYMBOL, FUNCTIONAL	68 212 05

**TABLE 6.6: RESUSCITAIRE® RADIANT WARMER CONTROLLER ASSEMBLY,
PARTS LIST**

(Sheet 2 of 3)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
29 AND 30		NOT USED	
31		POWER SUPPLY, SWITCHING, DC, 40W	81 300 15
32		NOT USED	
33		CABLE ASSY, DATA PORT	81 300 40
34		AC WIRE HARNESS ASSY	81 300 42
35		CABLE ASSY, D.C. POWER	81 300 44
36		CABLE ASSY, SKIN PROBE	81 300 46
37		CABLE ASSY, AUXILIARY SWITCH	81 300 48
38		CABLE ASSY, KEY PAD LOCK	81 300 50
39		WIRE ASSY, GROUND 2	81 300 52
40		CABLE ASSY, AC RCPT TO CKT BRKR	81 300 54
41 AND 42		NOT USED	
43		CONN ASSY, AC RCPT, MALE, 10A IEC	81 300 60
44		FAN ASSY	81 300 62
45 AND 46		NOT USED	
47		LABEL, DATA TAG	81 001 01
48 AND 49		NOT USED	
50		SCR, 4-40x1/4, TR PH SS	99 010 56
51		SCR, 4-40x7/16, TR PH SS	99 011 35
52		NOT USED	
53		SCR, 6-32X1/4, TR PH SS	99 022 72
54		SCR, 6-32X3/8, TR PH SS	99 023 26
55		NOT USED	
56		NUT, HEX, 6-32 CONL WSHR KEPS S	99 105 70
57		NOT USED	
58		WSHR, LK,INT, #4 SS	99 121 35
59		WSHR, LK EXT, #6 SS	99 122 20
60		WSHR, RIB SPRING S3 S ZI	99 999 72
61		NOT USED	
62		STANDOFF, HEX, 6-32 M/F 4.13L SS	99 116 87
63		STANDOFF, HEX, 6-32 M/F 4.785L SS	99 116 95
64		STANDOFF, HNG 6-32 ID .25OD 1.00L SS	99 116 89

**TABLE 6.6: RESUSCITAIRE® RADIANT WARMER CONTROLLER ASSEMBLY,
PARTS LIST**

(Sheet 3 of 3)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
65		STANDOFF, 6-32ID, .25HX.65L NYL	99 116 90
66		NOT USED	
67		LOCTITE, SCREWLOCK NO. 222	99 901 38
68 AND 69		NOT USED	
70		CONN, BD SPCR, ML DBL ROW 42 PSN	17 BS4 62
71		PROBE, AMBIENT TEMP	81 300 69
72		RING, RTNG, NAT 6100-37-ST-ZD	99 182 06
73		SUPPRESSOR, FERRITE 120 OHMS	17 585 58
74		WIRE TIE (4.50 LG)	12 995 00

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

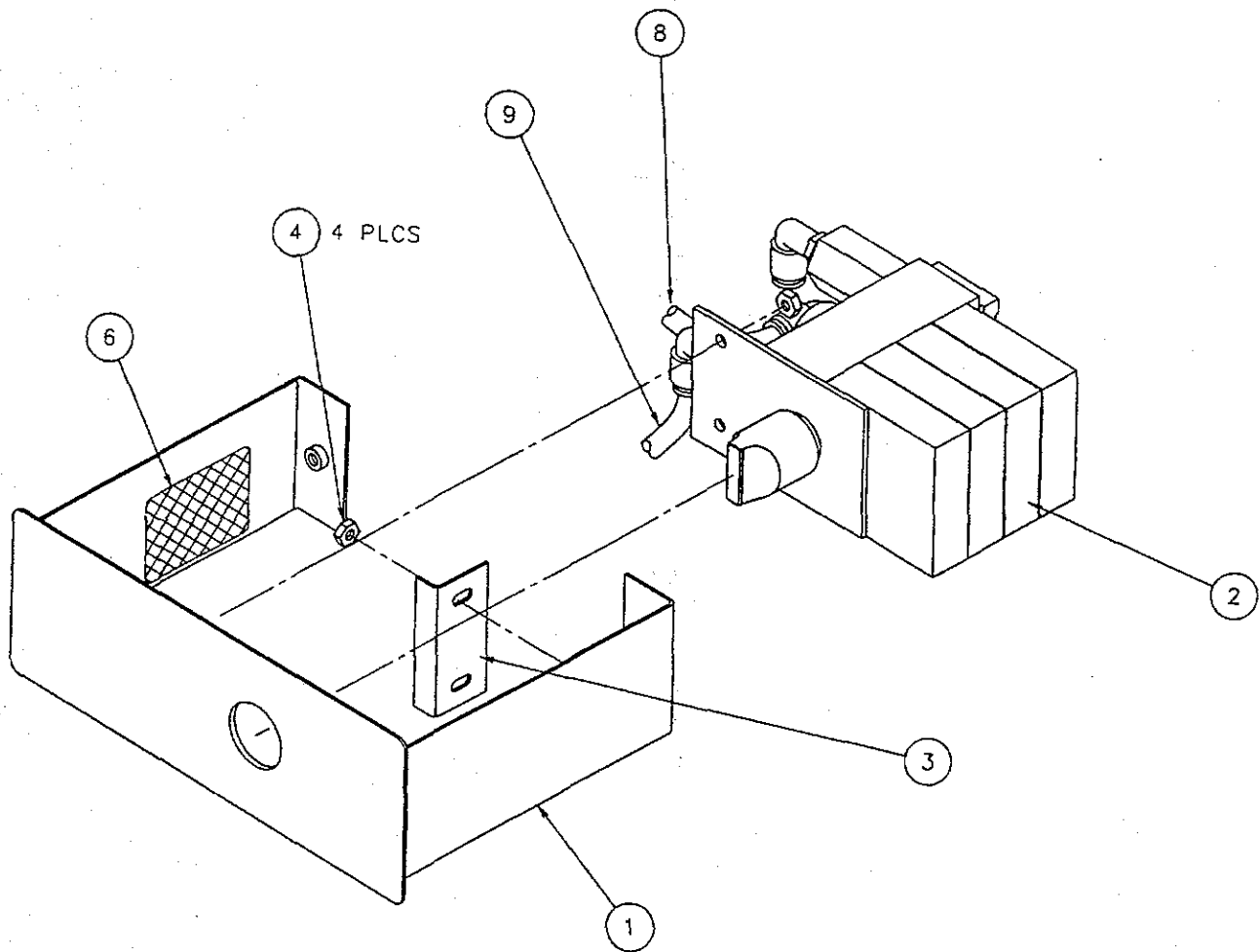


FIGURE 6.7: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
BLENDER ASSEMBLY
(SHEET 1 OF 1)

**TABLE 6.7: RESUSCITAIRE® RADIANT WARMER BLENDER ASSEMBLY,
PARTS LIST**

(Sheet 1 of 1)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		BLENDER ASSEMBLY	81 600 80
1		CHASSIS, MICROBLENDER  (Order label item 13. Refer to Figure 6.1.)	81 600 00
2		MICROBLENDER, LOW FLOW	81 600 40
3		CLIP, SUPPORT, MICROBLENDER	81 600 06
4		NUT, HEX, 8-32 KEPS S ZI	99 106 32
5		NOT USED	
6		LABEL, DATA TAG	81 000 38
7		NOT USED	
8		TUBING, POLY, 5/160OD .19ID 11.0L	81 501 38
9		TUBING, POLY, 5/160OD .19ID 14.0L	81 501 39

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		BLENDER BLANK ASSEMBLY	81 600 70
1		CHASSIS, MICROBLENDER  (Order label item 13. Refer to Figure 6.1.)	81 600 00
2 THRU 9		NOT USED	
10	NOT SHOWN	TUBING, POLY, 5/160OD .19ID 4.50L	81 501 36

**RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST**

Diagram to be supplied.

**FIGURE 6.8: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
RESUSCITATION MODULE (W/O THE AUTOBREATH INFANT RESUSCITATOR) ASSEMBLY
(SHEET 1 OF 1)**

**TABLE 6.8: RESUSCITAIRE® RADIANT WARMER
RESUSCITATION MODULE ASSEMBLY
(WITHOUT THE AUTOBREATH INFANT RESUSCITATOR),
PARTS LIST**

(Sheet 1 of 1).

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		RESUSCITATION MODULE ASSEMBLY, W/O AUTOBREATH (PRESSURE GAUGE PORT) (Order label item 14. Refer to Figure 6.1.)	81 400 72
		RESUSCITATION MODULE ASSEMBLY (NEW), W/O AUTOBREATH (PRESSURE GAUGE PORT) (Order label item 14. Refer to Figure 6.1.)	81 400 73

**RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST**

Diagram to be supplied.

**FIGURE 6.9: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
RESUSCITATION MODULE (WITH THE AUTOBREATH INFANT RESUSCITATOR) ASSEMBLY
(SHEET 1 OF 1)**

**TABLE 6.9: RESUSCITAIRE® RADIANT WARMER
RESUSCITATION MODULE ASSEMBLY
(WITH THE AUTOBREATH INFANT RESUSCITATOR),
PARTS LIST**

(Sheet 1 of 1)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		RESUSCITATION MODULE ASSEMBLY, WITH THE AUTOBREATH INFANT RESUSCITATOR (Order label item 14. Refer to Figure 6.1.)	81 400 80
		RESUSCITATION MODULE ASSEMBLY (NEW), WITH THE AUTOBREATH INFANT RESUSCITATOR (Order label item 14. Refer to Figure 6.1.)	81 400 81

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

NOTES:

1. CHAIN FOR ITEM 30 TO BE INSERTED THROUGH ITEM 1 BEFORE INSTALLING CAP.
2. COPPER TUBING, LETTER "A" IS PRE-FORMED.
3. COPPER TUBING, LETTER "B" IS FREE FORMED.
4. ITEM 11 IS TO BE ADJUSTED TO FIT FIRMLY AGAINST THE INSIDE OF THE UPPER POST AT FINAL ASSEMBLY.
5. ITEM 46 TO BE APPLIED TO ALL PIPE FITTINGS.
6. ITEM 47 TO BE APPLIED TO ALL THREADED FASTENERS.
7. FITTINGS MAY BE ROTATED FOR EASE OF TUBING CONNECTIONS.
8. POSITION IT.23 SO FLAT IS PERPENDICULAR TO SET SCREW.

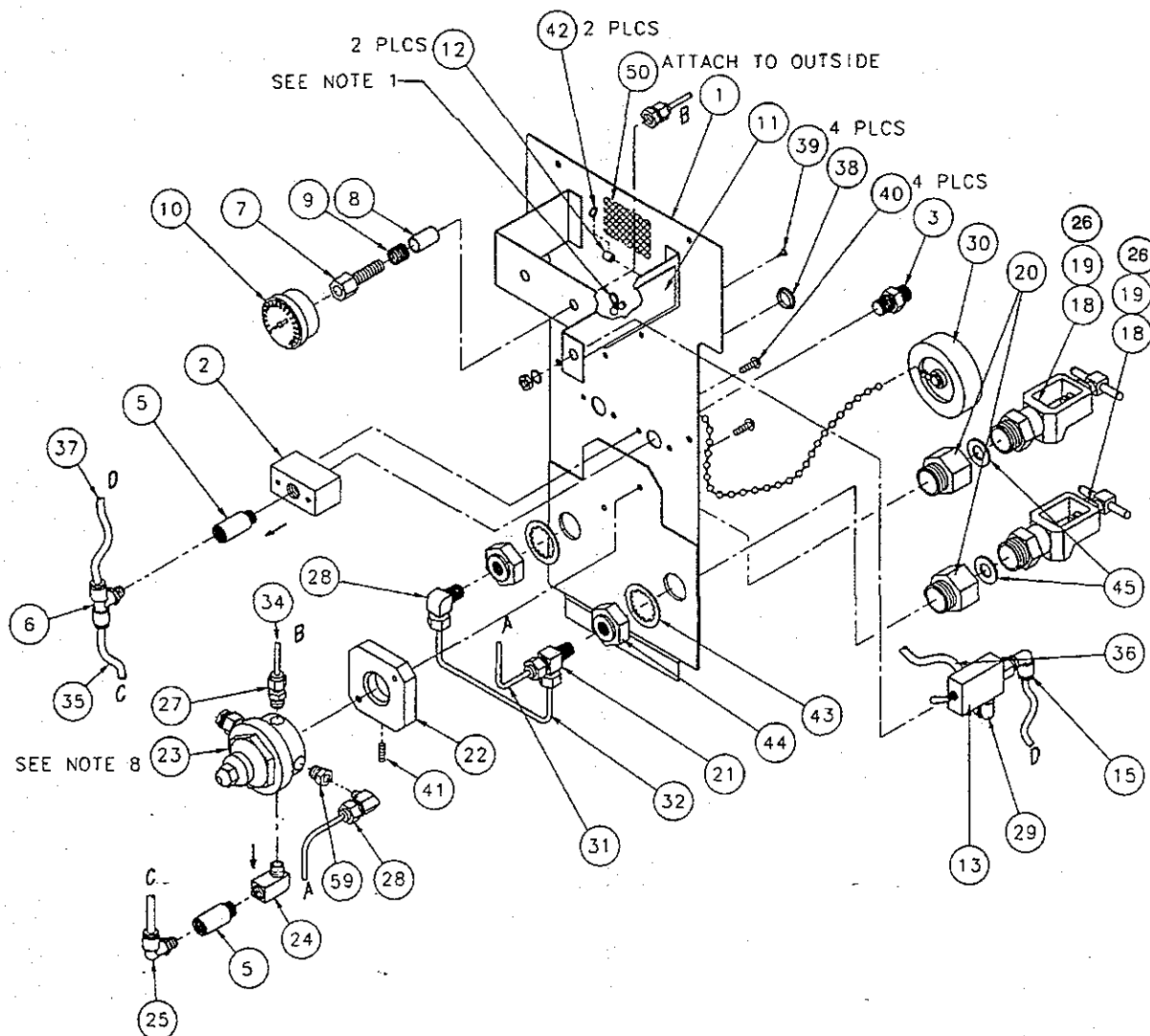


FIGURE 6.10: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
O₂ PIPELINE (DISS/NIST) WITH O₂ CYLINDER (PIN-INDEX) GAS SUPPLY ASSEMBLY
(SHEET 1 OF 1)

**TABLE 6.10: RESUSCITAIRE® RADIANT WARMER
O₂ PIPELINE (DISS/NIST) WITH O₂ CYLINDER (PIN-INDEX)
GAS SUPPLY ASSEMBLY, PARTS LIST**

(Sheet 1 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		GAS SUPPLY, O ₂ PIPELINE (DISS) O ₂ WITH CYLINDER (PIN-INDEX)	81 500 80
		GAS SUPPLY, O ₂ PIPELINE (NIST) O ₂ WITH CYLINDER (PIN-INDEX)	81 500 81
1		CHASSIS, GAS SUPPLY, W/CYL HOLES  (Order label items 2, 3, 9 and 15. Refer to Figure 6.1.)	81 500 01
2		MANIFOLD, PIPELINE	81 500 06
3		CONNECT, MALE, 1/4 MNPT TO 1/4 MNPT	68 239 27
		ADAPTER, 1/4NPT M-O2 NIST FEM	68 239 36
4		NOT USED	
5		CHECK VALVE, 1/4 MNPT TO 1/4 MNPT	81 500 43
6		FITTING, BR TEE, .250DTBG1/4MNPT	81 500 21
7		BULKHEAD, SWAGE, 125OD 1/8NPT F	81 500 32
8		SPACER, .315ID .375OD .480LG	81 500 03
9		SPRING, CPRS, .480OD, .045W, .875L	81 500 08
10		GAUGE, PRESSURE, 0-4000 PSI	81 500 36
11		BRACKET, PNEUMATIC VALVE	81 500 02
12		STANDOFF, 6-32ID, .25OD .25LG BR	99 116 91
13		VALVE, CONTROL, MINI 3-WAY	81 500 40
14		NOT USED	
15		ELBOW, .250OD TBG 1/8 NPT M	81 500 22
16 AND 17		NOT USED	
18		YOKE, OXYGEN (FOR CHECK VALVE)	81 500 46
19		CHECK VALVE	78 422 41
20		ADAPTER, YOKE, MALE TO FEMALE	78 440 15
21		TEE, RUN, SWAGE, 19OD TBG 1/8 NPT FM	81 501 72
22		BLOCK, MTG, PRESS. REGULATOR	81 500.05
23		REGULATOR, 65 PSI W/PRESS RELIEF	81 501 69
24		ELBOW, STREET, 90 DEG, 1/4NPT M/F	68 238 16
25		ELBOW, .250OD, 1/4 NPT M	81 500 23
26		NOT USED	

Change 1

**TABLE 6.10: RESUSCITAIRE® RADIANT WARMER
O₂ PIPELINE (DISS/NIST) WITH O₂ CYLINDER (PIN-INDEX)
GAS SUPPLY ASSEMBLY, PARTS LIST**

(Sheet 2 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
27		CONNECTOR, SWAGE, .125OD 1/4 NPT M	81 500 27
28		ELBOW, SWAGE, .190OD TBG 1/8 NPT M	81 500 30
29		ELBOW, .312OD TBG, 1/8MNPT	81 500 24
30		HANDLE, VALVE, REMOVABLE W/CHAIN	75 407 38
31		TUBING, FORMED, 3/16OD	81 501 08
32		TUBING, FORMED, 3/16OD	81 501 05
33		NOT USED	
34		TUBING, COPPER, .125OD .03WL 12.00LG	81 501 33
35		TUBING, PE, .25OD .170ID 7.00LG	81 501 11
36		TUBING, PE, .187ID .312OD 11.00 L	81 501 38
37		TUBING, PE, .25OD .170ID 7.75 LG	81 501 14
38		PLUG, BUTTON, NYLON, 0.62 HOLE	81 500 42
39		PLUG PNL FSTNG, .182 - .192, CANOE	68 602 02
40		SCR, 8-32x1/2, TR PH SS NYLOK	99 031 99
41		SCR, 10-32x1/2 SE SK SS CP	99 042 00
42		RING, RTNG, PUSH ON (5555-25)	99 181 48
43		WSHR, LK, INT, 7/8 S CA	99 129 85
44		NUT, HEX, JAM, 7/8-14 S ZI	99 114 50
45		WSHR, FL, .234ID x .81OD x .06T NYLN	78 423 30
46		TAPE, PIPE JNT SEALER, TEFLON, .50"	99 900 26
47		LOCTITE, SCREWLOCK NO. 222	99 901 38
48 AND 49		NOT USED	
50		LABEL, DATA TAG	81 000 38
51 THRU 58		NOT USED	
59		BUSHING, 1/4NPT-1/8NPT. BRASS	68 504 09

Change 1

This page is intentionally left blank.

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

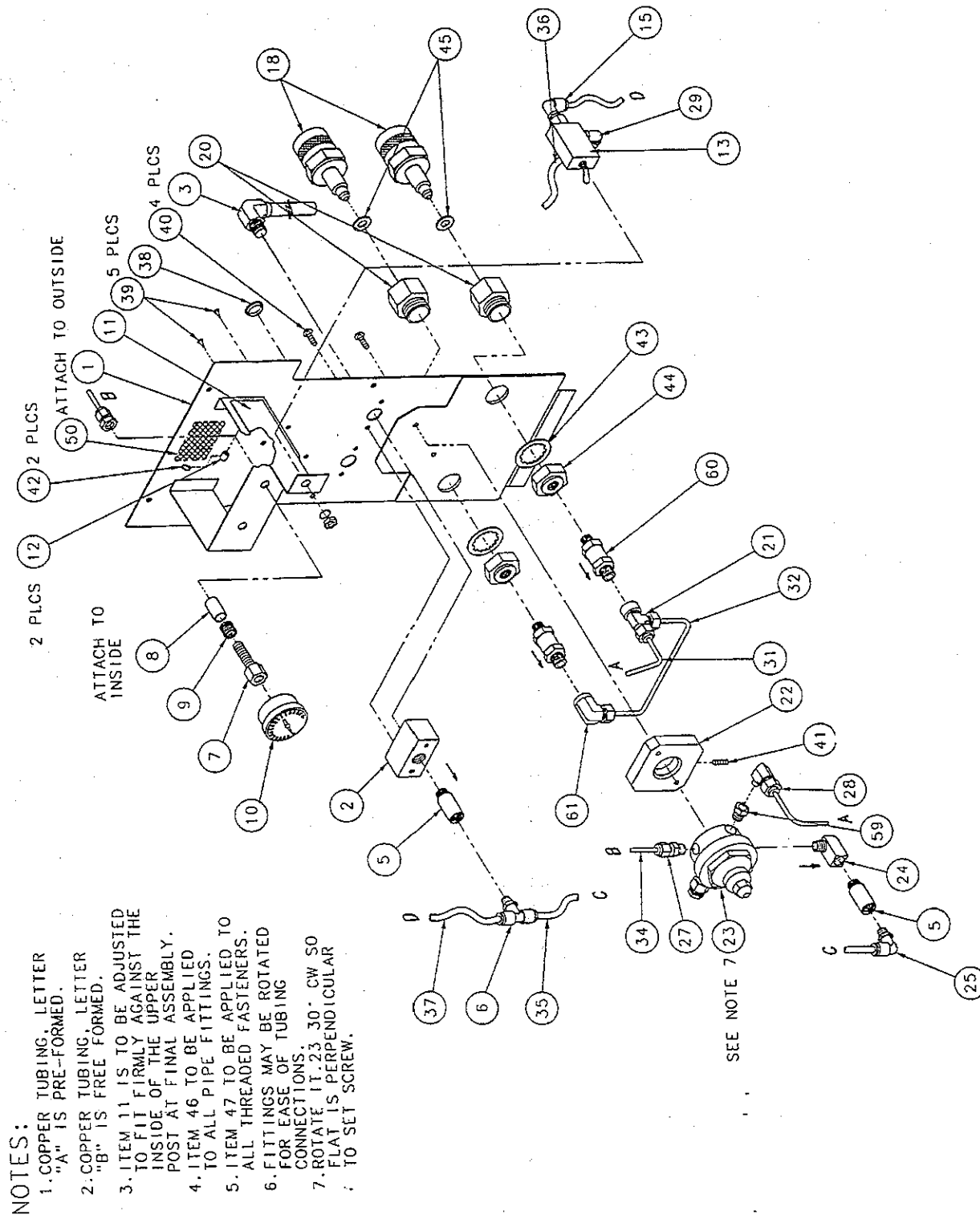


FIGURE 6.11: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
O₂ (WHITE) HOSE ASSEMBLY WITH O₂ CYLINDER (DIN FITTINGS) GAS SUPPLY ASSEMBLY,
(SHEET 1 OF 1)

Change 1

**TABLE 6.11: RESUSCITAIRE® RADIANT WARMER
O₂ (WHITE) HOSE ASSEMBLY WITH O₂ CYLINDER (DIN FITTINGS)
GAS SUPPLY ASSEMBLY, PARTS LIST**

(Sheet 1 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		GAS SUPPLY, O ₂ PIPELINE (WHITE HOSE), O ₂ CYLINDER (DIN FITTINGS)	81 500 83
1		CHASSIS, GAS SUPPLY, W/CYL HOLES  (Order label items 2, 3, 9 and 15. Refer to Figure 6.1.)	81 500 01
2		MANIFOLD, PIPELINE (WALL)	81 500 06
3		HOSE ASSY, OXYGEN, WHITE, 10 FT	81 501 53
4		NOT USED	
5		CHECK VALVE, 1/4NPT M IN, F OUT	81 500 43
6		TEE, BRANCH, .25OD TBG 1/4NPT M	81 500 21
7		BULKHEAD, SWAGE, 125OD 1/8NPT F	81 500 32
8		SPACER, .315ID .375OD .480LG	81 500 03
9		SPRING, CPRS, .480OD .045W .875L	81 500 08
10		GAUGE, PRESSURE, 0-4000 PSI	81 500 36
11		BRACKET, PNEUMATIC VALVE	81 500 02
12		STANDOFF, 6-32ID, .25OD, .25LG BR	99 116 91
13		VALVE, CONTROL, MINI 3-WAY	81 500 40
14		NOT USED	
15		FITTING, ELB, .250TBG 1/8MNPT	81 500 22
16 AND 17		NOT USED	
18		CONNECTOR, DIN, OXY CYL	81 501 22
19		NOT USED	
20		ADAPTER, HIGH PRESSURE CYLINDER	81 500 34
21		TEE, 19TBG SWG-1/8NPTF, 0.19TBG	81 501 72
22		BLOCK, MTG, PRESS. REGULATOR	81 500 05
23		REGULATOR, 65PSI W/PRESS RELIEF	81 500 69
24		ELBOW, STREET, 90 DEG, 1/4NPT M/F BR	81 500 37
25		FITTING, ELBOW, .250OD TBG 1/4NPT M	81 500 23

Change 2

**TABLE 6.11: RESUSCITAIRE® RADIANT WARMER
O₂ (WHITE) HOSE ASSEMBLY WITH O₂ CYLINDER (DIN FITTINGS)
GAS SUPPLY ASSEMBLY, PARTS LIST**

(Sheet 2 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
26		NOT USED	
27		CONNECTOR, SWAGE, .125OD 1/4NPT M	81 500 27
28		ELBOW, STREET, 90 DEG, 1/4 NPT M/F	81 500 30
29		ELBOW, .312OD TBG, 1/8NPT M	81 500 24
30		NOT USED	
31		TUBING, FORMED, 3/16OD	81 501 08
32		TUBING, FORMED, 3/16OD	81 501 05
33		NOT USED	
34		TUBING, COPPER, .125OD .03WL 12.0LG	81 501 33
35		TUBING, PE, .25OD .170ID 7.00LG	81 501 11
36		TUBING, PE, .187OD .312OD 11.00L	81 501 38
37		TUBING, PE, .25OD .170ID 7.75LG	81 501 14
38		PLUG BUTTON, NYLON, 0.62 HOLE	81 500 42
39		PLUG, PNL FSTNG, .182-.192, CANOE	68 602 02
40		SCR, 8-32X1/2, TR PH SS NYLOK	99 031 99
41		SCR, 10-32X1/2 SE SK SS CP	99 042 00
42		RING, RTNG,PUSH ON, SS (5555-25)	99 181 48
43		WSHR, LK, INT, 7/8 S CA	99 129 85
44		NUT, HEX, JAM, 7/8-14 S ZI	99 114 50
45		WSHR, FL, .190ID .562OD .062T NL	99 123 70
46		TAPE, PIPE JNT SEALER, TEFLON, .50"	99 900 26
47		LOCTITE, SCREWLOCK NO. 222	99 901 38
48 and 49		NOT USED	
50		LABEL, DATA TAG	81 000 38
51 thru 58		NOT USED	
59		BUSHING, 1/4NPT-1/8NPT, BRASS	68 504 09
60		CHECK VALVE, 1/4NPT-1/8NPT, BRASS	68 504 09
61		ELBOW, 90 DEG 0.19TBG-1/8NPT F BR	81 501 71

Change 1

This page is intentionally left blank.

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

NOTES:

1. CHAIN FOR ITEM 30 TO BE INSERTED THROUGH ITEM 1 BEFORE INSTALLING CAP.
2. COPPER TUBING, LETTERS A & D ARE PRE-FORMED.
3. COPPER TUBING, LETTERS B & C ARE FREE FORMED.
4. ITEM 11 IS TO BE ADJUSTED TO FIT FIRMLY AGAINST THE INSIDE OF THE UPPER POST AT FINAL ASSEMBLY.
5. ITEM 46 TO BE APPLIED TO ALL PIPE FITTINGS.
6. ITEM 47 TO BE APPLIED TO ALL THREADED FASTENERS.
7. FITTINGS MAY BE ROTATED FOR EASE OF TUBING CONNECTIONS.
8. POSITION IT. 23 WITH FLAT PERPENDICULAR TO SET SCREW IT. 41

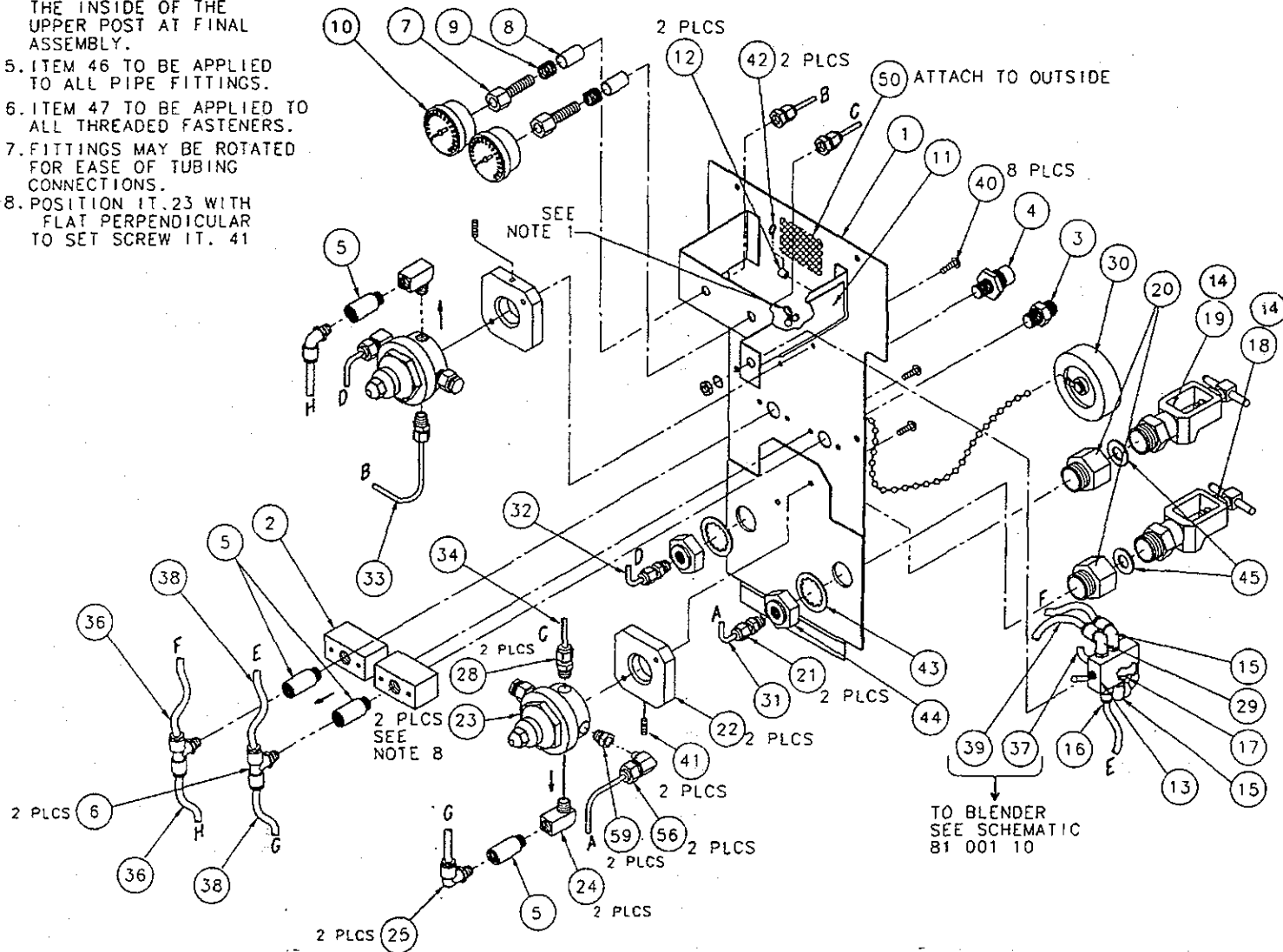


**FIGURE 6.12: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
O₂/AIR PIPELINE (DISS/NIST) WITH RESERVE (YOKE FITTING) GAS SUPPLY ASSEMBLY,
(SHEET 1 OF 1)**

**TABLE 6.12: RESUSCITAIRE® RADIANT WARMER
O₂/AIR PIPELINE (DISS/NIST) WITH O₂/AIR CYLINDER (PIN-INDEX)
GAS SUPPLY ASSEMBLY,
PARTS LIST**

(Sheet 1 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		GAS SUPPLY, O ₂ /AIR PIPELINE (DISS) WITH CYLINDER (PIN-INDEX)	81 500 85
		GAS SUPPLY, O ₂ /AIR PIPELINE (NIST) WITH CYLINDER (PIN-INDEX)	81 500 86
1		CHASSIS, GAS SUPPLY, W/CYL HOLES  (Order label items 2, 3, 9 and 15. Refer to Figure 6.1.)	81 500 01
2		MANIFOLD, PIPELINE	81 500 06
3		CONNECTOR, MALE, 1/4NPT-O ₂ DISS	68 239 27
		ADAPTER, 1/4NPT M-O ₂ NIST FEM	68 239 36
4		FITTING, DISS, MED AIR, 1/4MNPT	81 500 38
		ADAPTER, 1/4NPT M-AIR NIST FEM	81 500 47
5		CHECK VALVE, 1/4MNPT TO 1/4MNPT	81 500 43
6		FITTING, BR TEE, .25OD TBG 1/4MNPT	81 500 21
7		FITTING, BLKHD, SW.125OD, 1/8FNPT	81 500 32
8		SPACER, .315ID .375OD .480LG	81 500 03
9		SPRING, CMPR, .480OD, .045W, .875L	81 500 08
10		GAUGE, PRESSURE, 0-4000 PSI	81 500 36
11		BRACKET, PNEUMATIC VALVE	81 500 02
12		STANDOFF, 6-32ID, .25OD, .25LG BR	99 116 91
13		VALVE, CONTROL, MINI 5-WAY	81 500 41
14		NOT USED	
15		FITTING, ELB, .250TBG, 1/8MNPT	81 500 22
16		FITTING, CONN, .250OD 1/8MNPT	81 500 25
17		PIPE PLUG, 1/8NPT, HEX SOC HD, BR	78 456 08
18		YOKE, OXYGEN, MODIFIED	81 501 16
19		YOKE, MEDICAL AIR, MODIFIED	81 501 19
20		ADAPTER, YOKE, MALE TO FEMALE	78 440 15
21		FITTING, CONN, SW,.187OD, 1/8MNPT	81 500 28
22		BLOCK, MTG, PRESS. REGULATOR	81 500 05
23		REGULATOR, 65PSI W/PRESS RELIEF	81 501 69
24		ELBOW, STREET, 90 DEG, 1/4NPT M/F br	81 500 37
25		FITTING, ELBOW, .250OD, 1/4MNPT	81 500 23

Change 3

**TABLE 6.12: RESUSCITAIRE® RADIANT WARMER,
O₂/AIR PIPELINE (DISS/NIST) WITH O₂/AIR CYLINDER (PIN-INDEX)
GAS SUPPLY ASSEMBLY,
PARTS LIST**

(Sheet 2 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
26 AND 27		NOT USED	
28		FITTING, CONN, SW, .125OD, 1/4MNPT	81 500 27
29		FITTING, ELB, .312OD TBG, 1/8MNPT	81 500 24
30		HANDLE, VALVE, REMOVABLE W/CHAIN	75 407 38
31		TUBING, FORMED, 3/16OD	81 501 06
32		TUBING, FORMED, 3/16OD	81 501 07
33		TUBING, COPPER, .125OD .03WL 13.50LG	81 501 34
34		TUBING, COPPER, .125OD .03W 12.00LG	81 501 33
35		NOT USED	
36		TUBING, PE, .25OD .170ID 12.0LG	81 501 12
37		TUBING, PE, .25OD .170ID 21.0LG	81 501 13
38		TUBING, PE, .25OD .170ID 7.75LG	81 501 14
39		TUBING, POLY, 5/16OD .19ID 7.5L	81 501 37
40		SCR, 8-32X1/2, TR PH SS NYLOK	99 031 99
41		SCR, 10-32X1/2, SE SK SS CP	99 042 00
42		RING, RTNG, PUSH ON (5555-25)	99 181 48
43		WSHR, LK, INT, 7/8 S CA	99 129 85
44		NUT, HEX, JAM, 7/8-14 S ZI	99 114 50
45		WSHR, FL, .234ID x .810D x .06T NYLON	78 423 30
46		TAPE, PIPE JNT SEALER, TEFLON, .50"	99 900 26
47		LOCTITE, SCREWLOCK NO. 222	99 901 38
48 AND 49		NOT USED	
50		LABEL, DATA TAG	81 000 38
51 THRU 55		NOT USED	
56		ELBOW, SWAGE, .190D TBG 1/8NPT M	81 500 30
57 AND 58		NOT USED	
59		BUSHING, 1/4NPT-1/8NPT, BRASS	68 504 09

Change 1

This page is intentionally left blank.

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

NOTES:

1. COPPER TUBING, LETTERS A & D ARE PRE-FORMED.
2. COPPER TUBING, LETTERS B & C ARE FREE FORMED.
3. ITEM 11 IS TO BE ADJUSTED TO FIT FIRMLY AGAINST THE INSIDE OF THE UPPER POST AT FINAL ASSEMBLY.

4. ITEM 46 TO BE APPLIED TO ALL PIPE FITTINGS.
5. ITEM 47 TO BE APPLIED TO ALL THREADED FASTENERS.
6. FITTINGS MAY BE ROTATED FOR EASE OF TUBING CONNECTIONS.
7. POSITION IT.23 WITH FLAT PERPENDICULAR TO SET SCREW IT.41.

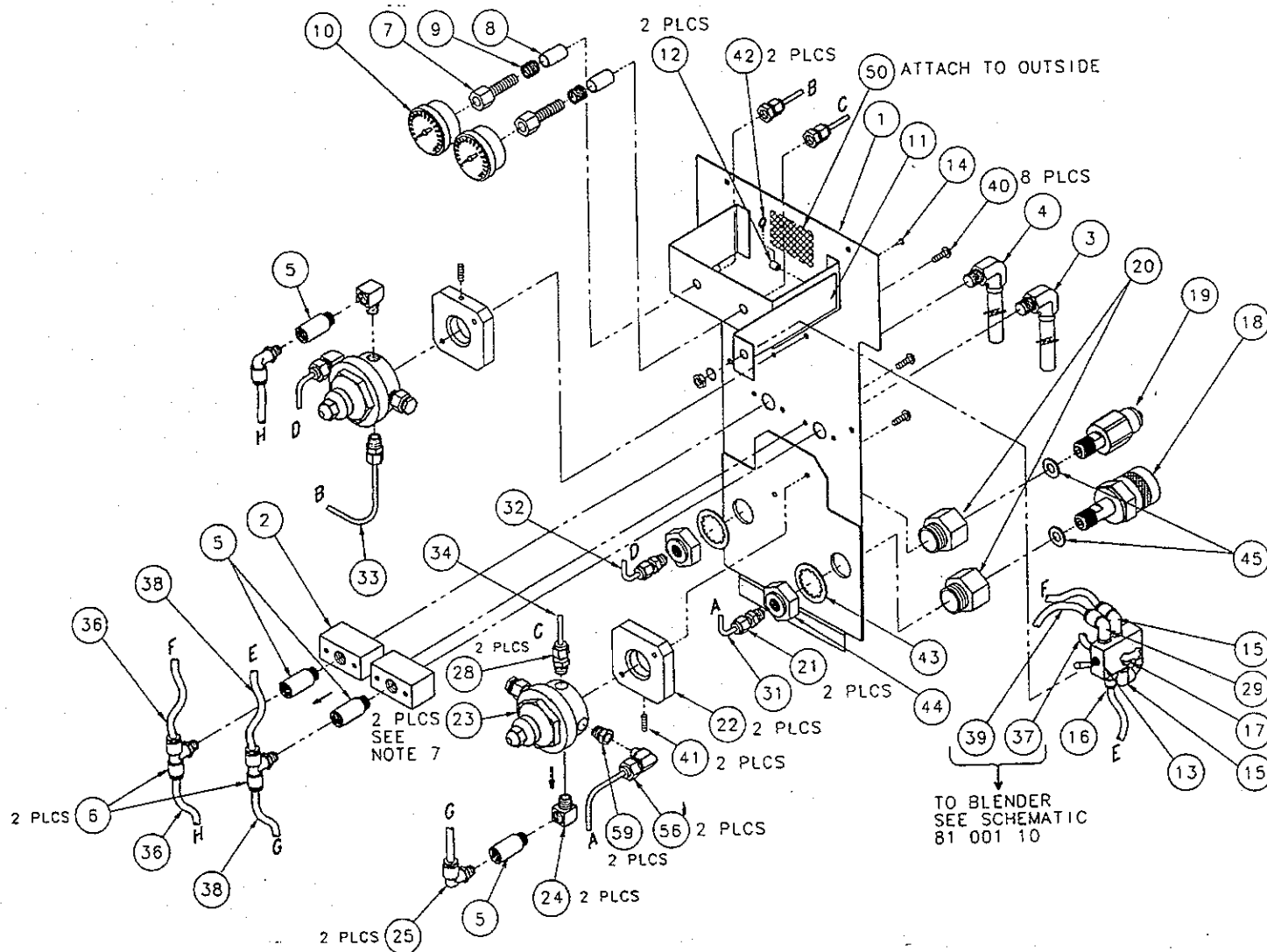


FIGURE 6.13: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
O₂/AIR HOSE ASSEMBLY WITH CYLINDER (DIN FITTINGS) GAS SUPPLY ASSEMBLY,
(SHEET 1 OF 1)

**TABLE 6.13: RESUSCITAIRE® RADIANT WARMER
O₂/AIR HOSE WITH CYLINDER (DIN FITTINGS) GAS SUPPLY ASSEMBLY,
PARTS LIST**

(Sheet 1 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
		GAS SUPPLY, O ₂ /AIR PIPELINE (BLUE/YELLOW HOSES), CYLINDER DIN FITTINGS	81 500 88
		GAS SUPPLY, O ₂ /AIR PIPELINE (GREEN/YELLOW HOSES), CYLINDER DIN FITTINGS	81 500 90
		GAS SUPPLY, O ₂ /AIR PIPELINE (BLUE/BLACK HOSES), CYLINDER DIN FITTINGS	81 501 91
1		CHASSIS, GAS SUPPLY, W/CYL HOLES  (Order label items 2, 3, 9 and 15. Refer to Figure 6.1.)	81 500 01
2		MANIFOLD, PIPELINE (WALL)	81 500 06
3		HOSE ASSY, OXYGEN, BLUE, 10 FT	81 501 52
		HOSE ASSY, OXYGEN, GREEN, 10 FT	81 501 55
4		HOSE ASSY, MED AIR, YELLOW, 10 FT	81 501 54
		HOSE ASSY, MED AIR, BLACK, 10 FT	81 501 51
5		CHECK VALVE, 1/4NPT F IN, M OUT	81 500 43
6		TEE, BRANCH, .25OD TBG 1/4NPT M	81 500 21
7		BULKHEAD, SWAGE, 125OD 1/8NPT F	81 500 32
8		SPACER, .315ID .375OD .480LG	81 500 03
9		SPRING, CPRS, .480OD .045W .875L	81 500 08
10		GAUGE, PRESSURE, 0-4000 PSI	81 500 36
11		BRACKET, PNEUMATIC VALVE	81 500 02
12		STANDOFF, 6-32ID, .25OD, .25LG BR	99 116 91
13		VALVE, CONTROL, MINI 5-WAY	81 500 41
14		PLUG, PNL FSTNG, .182 - .192, CANOE	68 602 02
15		FITTING, ELB, .250TBG 1/8MNPT	81 500 22
16		FITTING, CONN, .250OD 1/8MNPT	81 500 25
17		PIPE PLUG, 1/8NPT, HEX SOC HD, BR	78 456 08
18		CONNECTOR, DIN, OXY CYL	81 501 22
19		CONNECTOR, DIN, AIR CYL	81 501 23
20		ADAPTER, HIGH PRESSURE CYLINDER	81 500 34
21		FITTING, CONN, SW, .187OD, 1/8NPTM	81 500 28
22		BLOCK, MTG, PRESS. REGULATOR	81 500 05
23		REGULATOR, 65PSI W/PRESS RELIEF	81 501 69
24		ELBOW, STREET, 90 DEG, 1/4NPT M/F br	81 500 37
25		FITTING, ELBOW, .250OD TBG 1/4NPT M	81 500 23

Change 2

**TABLE 6.13: RESUSCITAIRE® RADIANT WARMER,
O₂/AIR HOSE WITH CYLINDER (DIN FITTINGS) GAS SUPPLY ASSEMBLY,
PARTS LIST**

(Sheet 2 of 2)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
26 and 27		NOT USED	
28		CONNECTOR, SWAGE, .125OD 1/4NPT M	81 500 27
29		ELBOW, .312OD TBG, 1/8NPT M	81 500 24
30		NOT USED	
31		TUBING, FORMED, 3/16OD	81 501 06
32		TUBING, FORMED, 3/16OD	81 501 07
33		TUBING, COPPER, .125OD .03WL 13.50LG	81 501 34
34		TUBING, COPPER, .125OD .03WL 12.0LG	81 501 33
35		NOT USED	
36		TUBING, POLY, .25OD .170ID 12.0LG	81 501 12
37		TUBING, POLY, .25OD .170ID 21.0LG	81 501 13
38		TUBING, POLY, .25OD .170ID 21.0LG	81 501 14
39		TUBING, POLY, .187ID .312OD 7.5LG	81 501 37
40		SCR, 8-32X1/2, TR PH SS NYLOK	99 031 99
41		SCR, 10-32X1/2 SE SK SS CP	99 042 00
42		RING, RTNG, PUSH ON(5555-25)	99 181 48
43		WSHR, LK, INT, 7/8 S CA	99 129 85
44		NUT, HEX, JAM, 7/8-14 S ZI	99 114 50
45		WSHR, FL, .190ID .562OD .062T NL	99 123 70
46		TAPE, PIPE JNT SEALER, TEFLON, .50"	99 900 26
47		LOCTITE, SCREWLOCK NO. 222	99 901 38
48 and 49		NOT USED	
50		LABEL, DATA TAG	81 000 38
51 thru 55		NOT USED	
56		ELBOW, SWAGE, .190OD TUBE, 1/8NPT M	81 500 30
57 and 58		NOT USED	
59		BUSHING, 1/4NPT-1/8NPT, BRASS	68 504 09

Change 1

This page is intentionally left blank.

RESUSCITAIRE® RADIANT WARMER

MODEL RW82-1

PARTS LIST

NOTES:

1. ASSEMBLY OF BOARD TO CONFORM WITH SPEC. NO. 00 945 12.
2. MARK ASSEMBLY NUMBER AND REVISION LEVEL IN AREA SHOWN WITH .09 ±.03 HIGH CHARACTERS USING A PERMANENT SMUDGE-PROOF INK, COLOR TO CONTRAST WITH BOARD MATERIAL.
3. COMPONENT LEAD AND SOLDER BUILD UP NOT TO EXCEED .06 MAX. FROM SURFACE OF BOARD.
4. CAUTION: ASSEMBLY CONTAINS STATIC SENSITIVE DEVICES, SPECIAL HANDLING REQUIRED.
5. * INDICATES VARIABLES AT ITEM NO'S. 11,12,27,28,29,31,34,60, 61,66,80.
6. SPACER UNDER 3 LED'S MUST BE PUSHED IN DIRECTION OF ARROW TO ASSURE THAT LED "DS7" HAS BOTTOM CONTACT OF BOTH ITS LEGS WITH SPACER; SPACER MUST BE HELD IN THIS POSITION DURING SOLDERING.

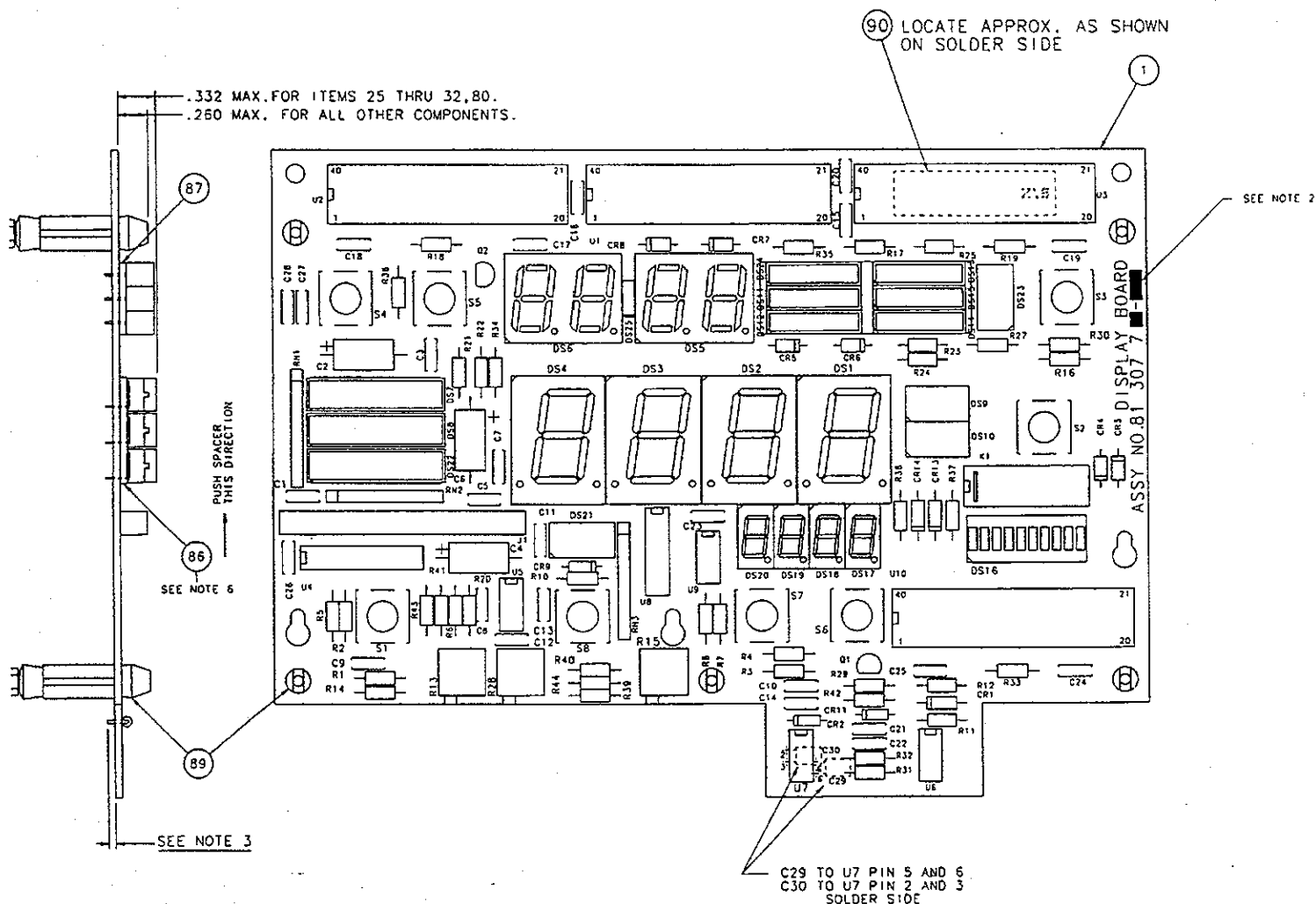


**FIGURE 6.14: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
PCB1 DISPLAY BOARD ASSEMBLY
(SHEET 1 OF 1)**

**TABLE 6.14: RESUSCITAIRE® RADIANT WARMER
PCB1 DISPLAY BOARD ASSEMBLY,
PARTS LIST**

(Sheet 1 of 1)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
-------------	---------------------	-------------	----------------

THE MAIN BOARD IS A MULTI-LAYERED BOARD. ATTEMPTING TO REPLACE ANY
COMPONENTS MAY RESULT IN DAMAGE TO THE BOARD.

PCB1 ASSY, DISPLAY, W/BABY MODE X81 307 75T

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

NOTES

1. ASSEMBLY OF BOARD TO CONFORM WITH SPEC. NO. 00 945 12.
2. MARK ASSEMBLY NUMBER AND REVISION LEVEL IN AREA SHOWN WITH .09 ±.03 HIGH CHARACTERS USING PERMANENT SMUDGE-PROOF INK, COLOR TO CONTRAST WITH BOARD MATERIAL.
3. COMPONENT LEAD AND SOLDER BUILD UP NOT TO EXCEED .06 MAX. FROM SURFACE OF BOARD.
4. CAUTION: ASSEMBLY CONTAINS STATIC SENSITIVE DEVICES, SPECIAL HANDLING REQUIRED.

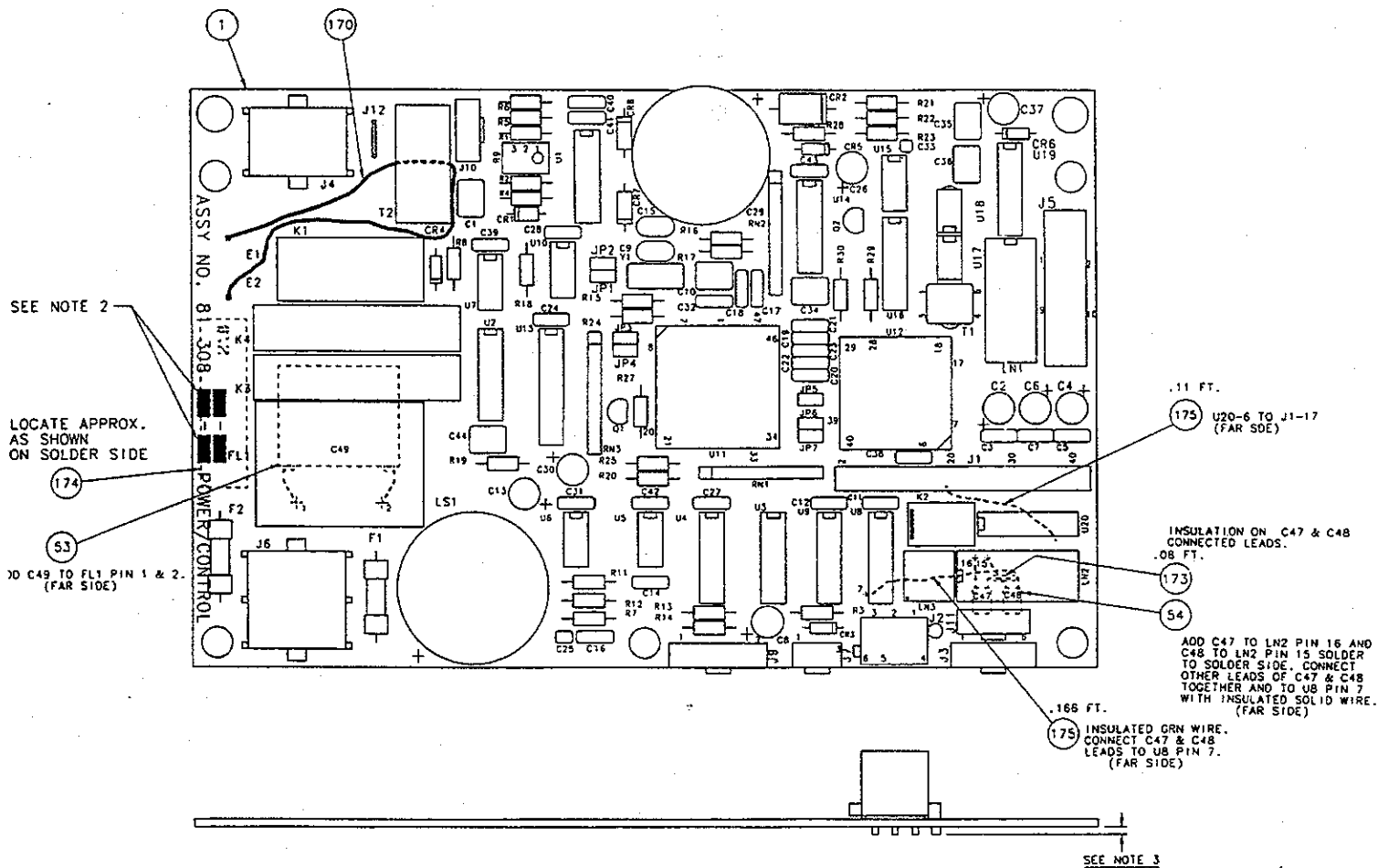


FIGURE 6.15: PARTS LOCATION DIAGRAM, RESUSCITAIRE® RADIANT WARMER
PCB2 POWER AND CONTROL BOARD ASSEMBLY
(SHEET 1 OF 1)

**TABLE 6.15: RESUSCITAIRE® RADIANT WARMER
PCB2 POWER AND CONTROL BOARD ASSEMBLY,
PARTS LIST**

(Sheet 1 of 1)

ITEM NO.	REFERENCE DESIG.	DESCRIPTION	PART NUMBER
-------------	---------------------	-------------	----------------

**THE MAIN BOARD IS A MULTI-LAYERED BOARD. ATTEMPTING TO REPLACE ANY
COMPONENTS MAY RESULT IN DAMAGE TO THE BOARD.**

PCB2 ASSY, POWER AND CONTROL BOARD, 120 V	X81 308 70T
PCB2 ASSY, POWER AND CONTROL BOARD, 220-240 V	X81 308 80T
PCB2 ASSY, POWER AND CONTROL BOARD, 220 V	X81 308 83T

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

This page is intentionally left blank.

SECTION 7 DIAGRAMS

7.1 GENERAL

This section provides the schematic diagrams for the Resuscitaire® Radiant Warmer.

This page is intentionally left blank.



FIGURE 7.1: RESUSCITAIRE® RADIANT WARMER, F

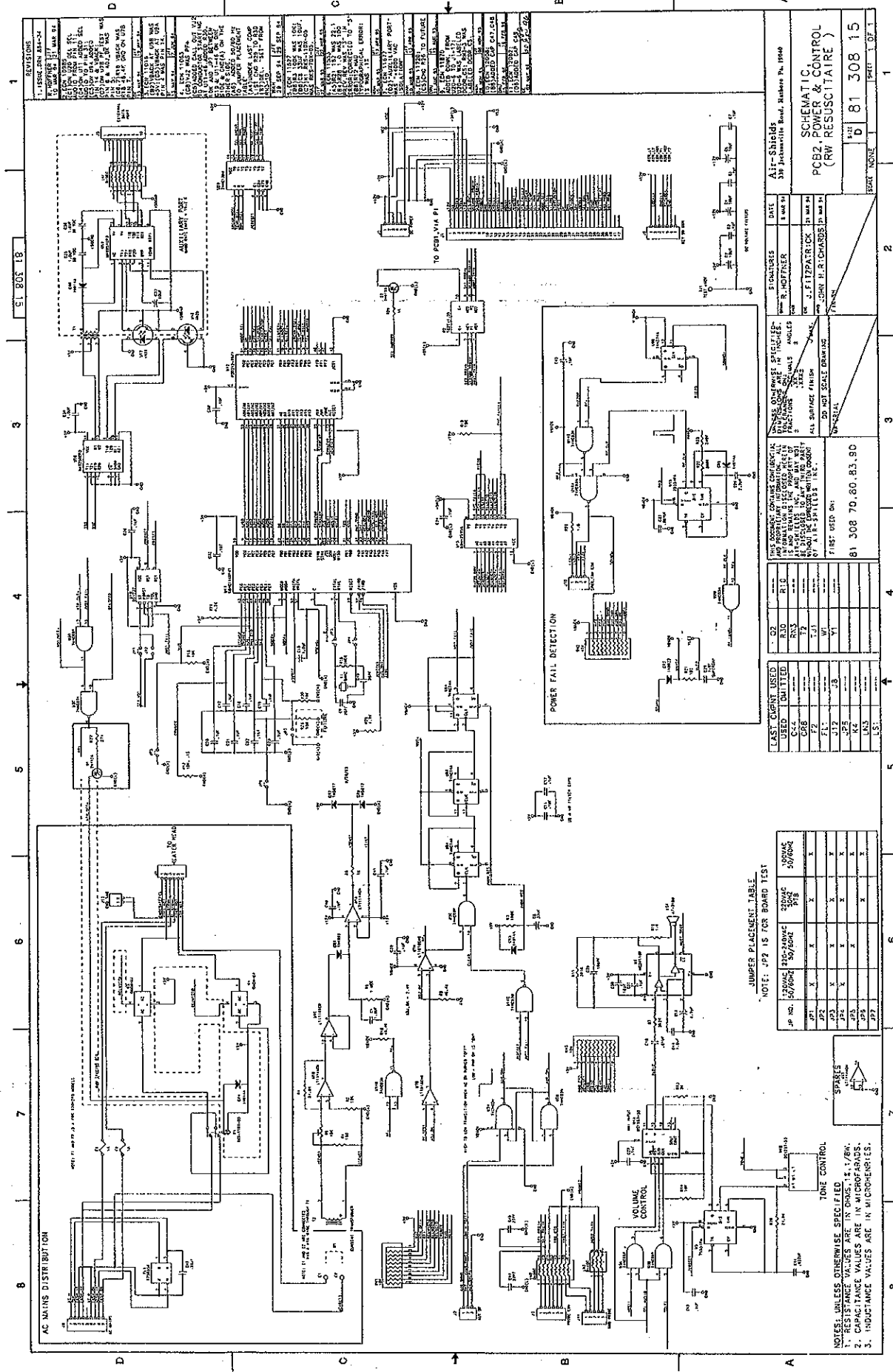


FIGURE 7-2: RESUSCITAIRE® RADIANT WARMER, P1


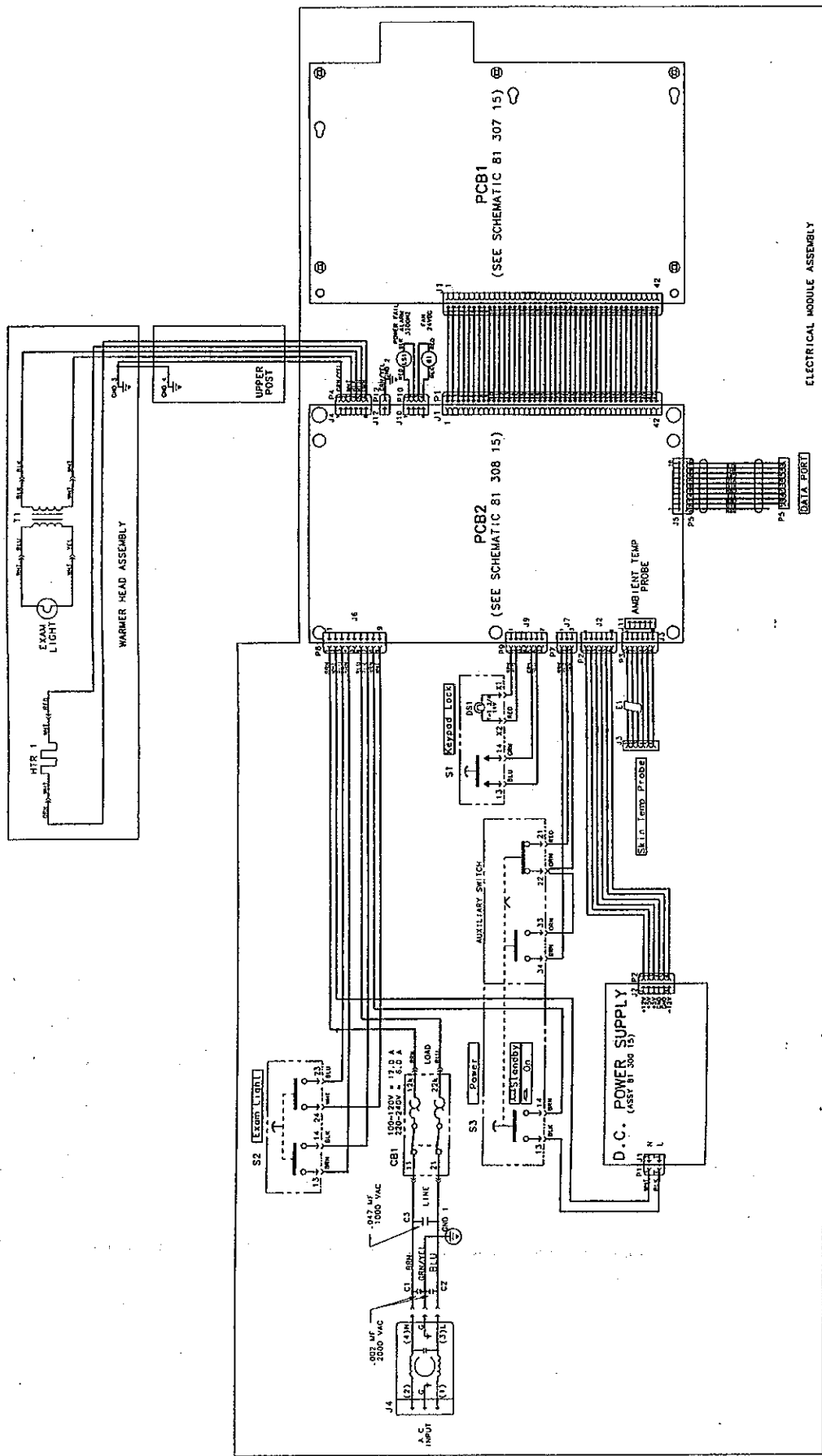



FIGURE 7.3 RESUSCITAIRE® RADIANT WARMER, INTERCONNECT, DIAGRAM SERIES 00
(SHEET 1 OF 1)
Change 1



FIGURE 7.3A RESUSCITAIRE® RADIANT WARMER, INTERCONNECT, DIAGRAM SERIES 01
(SHEET 1 OF 1)
Change 1

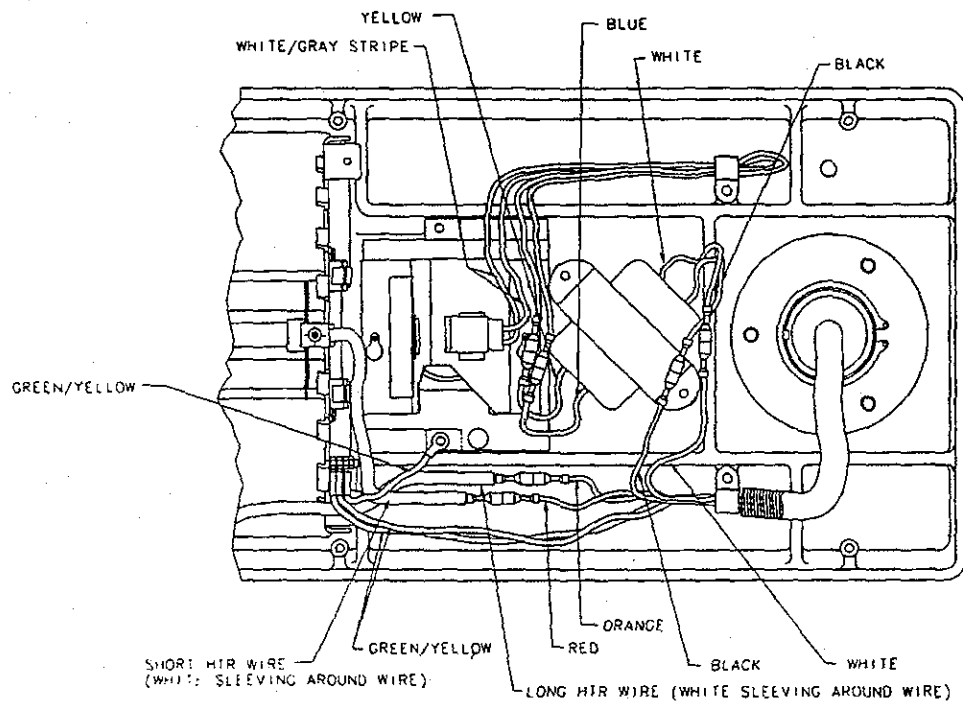


TOP VIEW
SHOWING
WIRE ROUTING THROUGH
CABLE CLAMPS
SEE SHEET 1 FOR ITEM NUMBERS

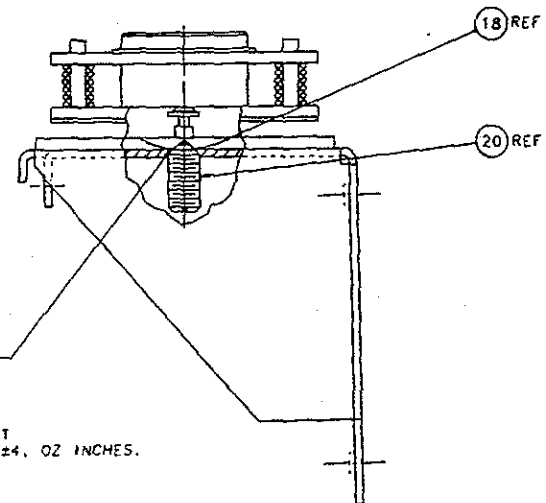

SIDE VIEW
VIEW INDICATES THE PROPER
INSERTION DEPTH
OF ITEM 20


FIGURE 7.4: RESUSCITAIRE® RADIANT WARMER, WARMER HEAD ASSEMBLY WIRING
SCHEMATIC SERIES 00
(SHEET 1 OF 1)
Change 1

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
DIAGRAMS

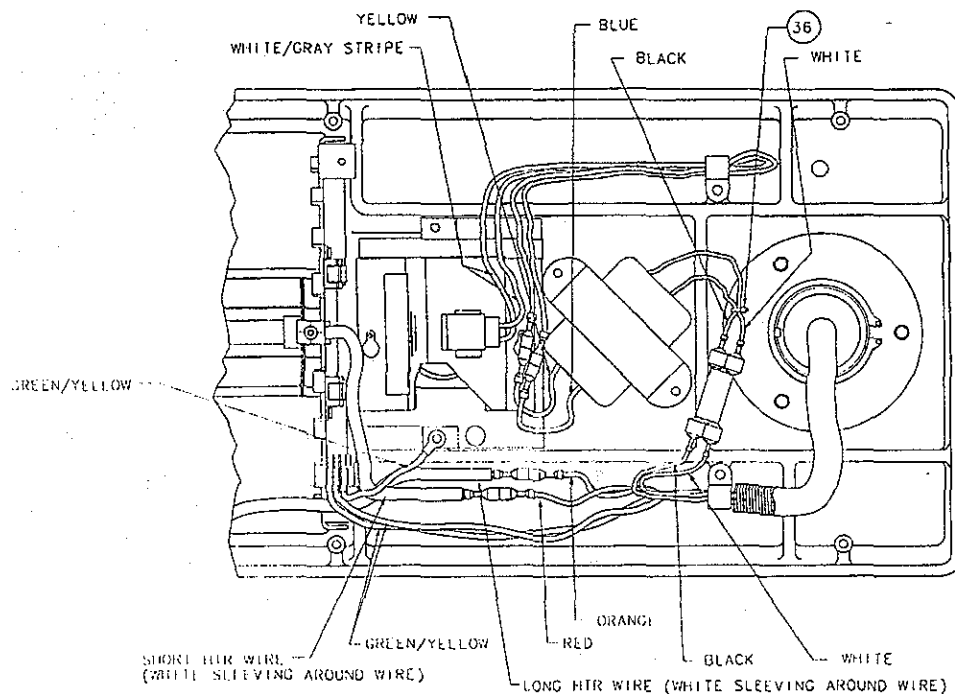


TOP VIEW
SHOWING
WIRE ROUTING THROUGH
CABLE CLAMPS

SEE SHEET 1 FOR ITEM NUMBERS

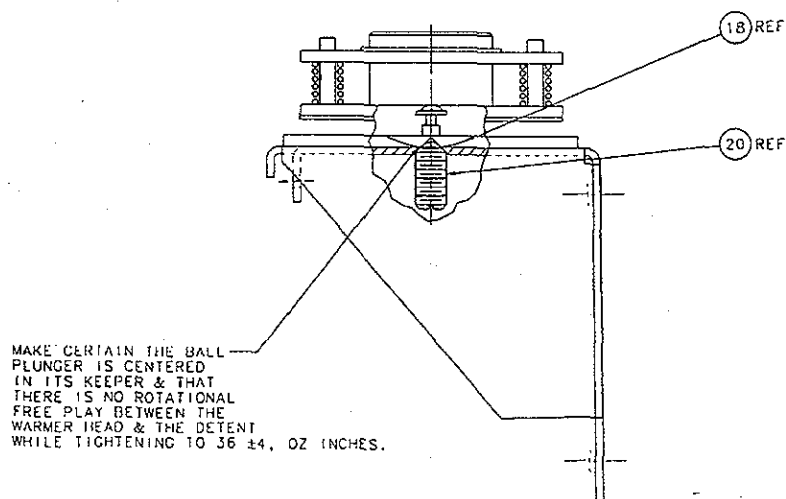

SIDE VIEW
VIEW INDICATES THE PROPER
INSERTION DEPTH
OF ITEM 20

FIGURE 7.4A: RESUSCITAIRE® RADIANT WARMER, WARMER HEAD ASSEMBLY WIRING
SCHEMATIC SERIES 01
(SHEET 1 OF 1)
Change 1

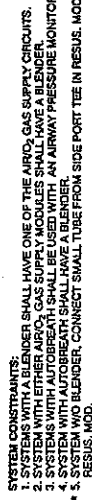

RECOMMENDED WITH AUTOBREATH

**FIGURE 7.5: RESUSCITAIRE® RADIANT WARMER
PNEUMATIC SYSTEM
(WITH THE AUTOBREATH INFANT RESUSCITATOR)
SCHEMATIC
(SHEET 1 OF 1)**

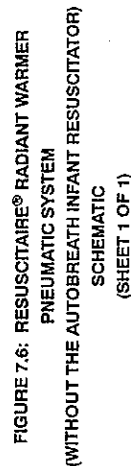


FIGURE 7.6: RESUSCITAIRE® RADIANT WARMER
PNEUMATIC SYSTEM
(WITHOUT THE AUTOBREATH INFANT RESUSCITATOR)
SCHEMATIC
(SHEET 1 OF 1)

RESUSCITAIRE® RADIANT WARMER
MODEL RW82-1
PARTS LIST

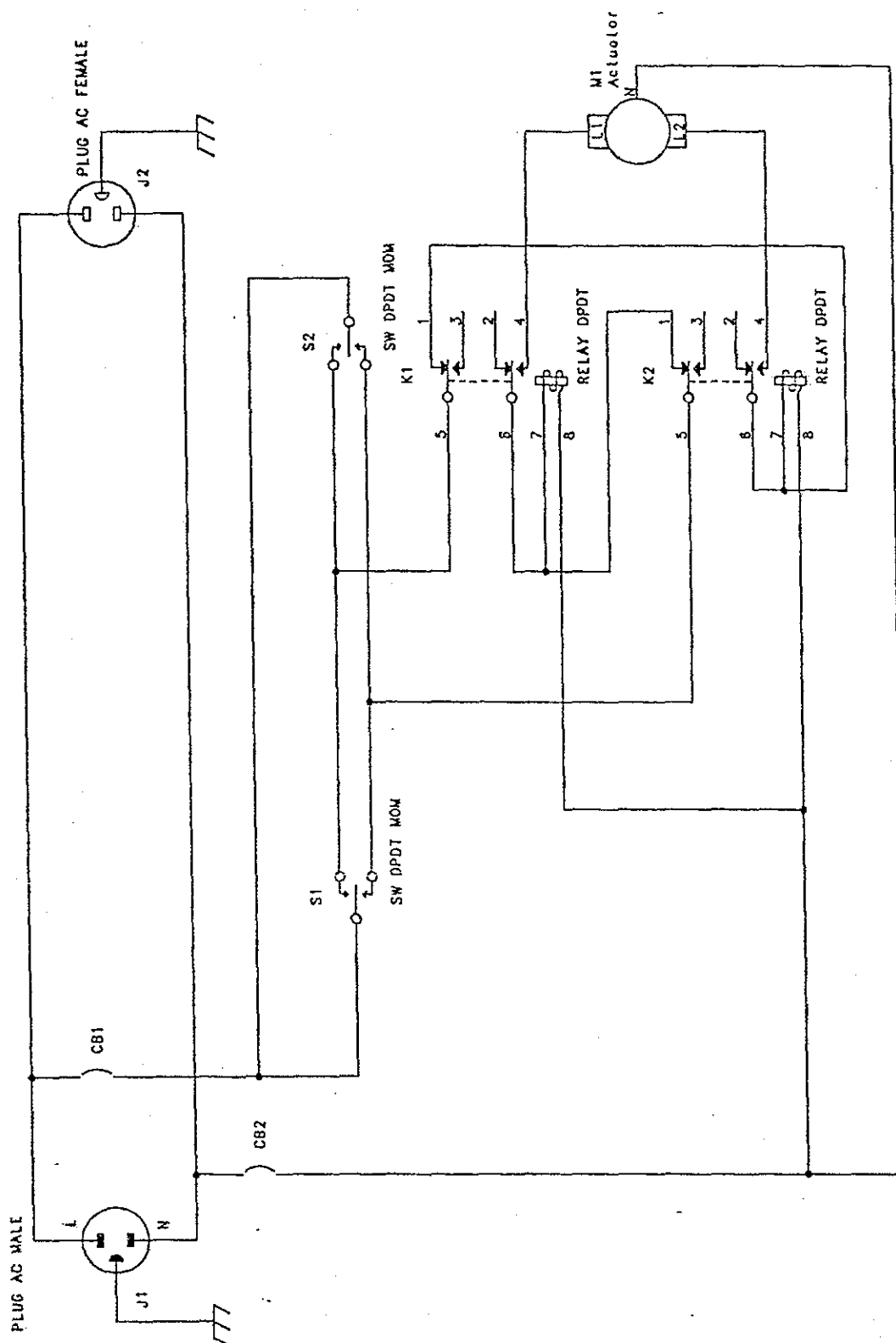


FIGURE 7.7 SCHEMATIC DIAGRAM, POWER ACTUATOR

Change 2

Hill-Rom Air-Shields®
A HILLENBRAND INDUSTRY

Hill-Rom Air-Shields reserves the right to make changes without notice in design, specifications, and models. The only warranty Hill-Rom Air-Shields makes is the express written warranty extended on the sale or rental of its products.